



Europäisches Patentamt
European Patent Office
Office européen des brevets



(11) Publication number:

0 431 149 B1

(12)

EUROPEAN PATENT SPECIFICATION

(45) Date of publication of patent specification: 06.09.95 (51) Int. Cl.⁸: C12Q 1/68

(21) Application number: 90913802.6

(22) Date of filing: 31.05.90

(86) International application number:
PCT/US90/03004

(87) International publication number:
WO 90/15157 (13.12.90 90/28)

(54) UNIVERSAL EUBACTERIA NUCLEIC ACID PROBES AND METHODS.

(30) Priority: 31.05.89 US 359158

(43) Date of publication of application:
12.06.91 Bulletin 91/24

(45) Publication of the grant of the patent:
06.09.95 Bulletin 95/36

(84) Designated Contracting States:
AT BE CH DE DK FR GB IT LI LU NL SE

(56) References cited:
EP-A- 245 129
EP-A- 250 662
EP-A- 277 237
EP-A- 0 272 009

SCIENCE, vol. 243, March 1989, E.F. DeLONG
et al., p. 1360&NUM;

DIALOG INFORMATIONAL SERVICES, file
154, medline 83-90, acc.no. 06930668, med-
line acc.no. 89232668; K. CHEN et al.&NUM;

(73) Proprietor: AMOCO CORPORATION
200 East Randolph Drive
P.O. Box 87703
Chicago
Illinois 60680-0703 (US)

(72) Inventor: LANE, David, J.
9 Oriole Drive
Milford, MA 01757 (US)
Inventor: SHAH, Jyotsna
13 Bates Drive
Nashua, NH 03060 (US)
Inventor: BUHARIN, Amelia
7 Pond Street
Framingham, MA 01701 (US)
Inventor: WEISBURG, William, G.
3 Jillson Circle
Milford, MA 01757 (US)

(74) Representative: Sheard, Andrew Gregory et al
Kilburn & Strode
30, John Street
London WC1N 2DD (GB)

Note: Within nine months from the publication of the mention of the grant of the European patent, any person may give notice to the European Patent Office of opposition to the European patent granted. Notice of opposition shall be filed in a written reasoned statement. It shall not be deemed to have been filed until the opposition fee has been paid (Art. 99(1) European patent convention).

DIALOG INFORMATIONAL SERVICES, file
154, medline 83-90, acc.no. 06470172, med-
line acc.no. 88115172; S.J. GIOVANNONI et
al.&NUM;

Description

This invention relates to detection of bacteria in clinical and other samples. Methods for the detection of bacteria in ordinarily aseptic bodily tissues or fluids such as blood, urine, and cerebrospinal fluid - where the presence of any bacterium may be life threatening are of particular importance. The present invention provides nucleic acid probes and compositions along with methods for their use for the specific detection of any bacterium in such samples.

EP-A-0272009 (Hogan *et al*, derived from WO-A-8802957) describes a number of nucleic acid probes, that hybridise to bacterial RNA but not to yeast or human rRNA, and a method for their preparation, that can be used in hybridisation assays for detecting non-viral organisms.

EP-A-0245129 (Institut Pasteur and INSERM) refers to oligonucleotide probes capable of hybridising to bacterial rRNA and their use in hybridisation assays for the detection of bacteria.

Background of the Invention

The term "eubacteria" as used herein, refers to the group of prokaryotic organisms (bacteria) as described in, for example, Bergey's Manual of Systematic Bacteriology (N.R. Krieg and J.G. Holt, ed., 1984, Williams & Wilkins, Baltimore). As a group, the eubacteria comprise all of the bacteria which are known to cause disease in humans or animals and are of most concern with respect to detection.

The only other described group of bacteria, the archaeobacteria, are biologically and genetically distinct from the eubacteria (C.R. Woese, Scientific American, 1981, Volume 244, pages 90-102). Archaeobacteria as a group occupy a variety of "extreme" environments such as hot springs, strongly oxygen-depleted muds, salt brines, etc., which generally do not support the growth of eubacteria. There are no known archaeobacterial pathogens and, consequently, their detection is of little clinical significance.

Eukaryotic organisms comprise the third fundamental genetic lineage which, together with the eubacteria and archaeobacteria, include all known life forms (Figure 1). Eukaryotes include humans, animals, plants and a host of organizationally less complex, free-living and parasitic "protists," including: protozoans, fungi, ciliates, etc. In a clinical context, it is particularly important that eubacteria be distinguished from eukaryotic, e.g. fungal and protozoan, infections which may present the same symptoms but require a significantly different regime of antimicrobial or chemo-therapy. These genetic distinctions thus are clinically significant from the point of view of diagnosis and antimicrobial chemotherapy.

It is an aspect of the present invention to provide nucleic acid probes which discriminate between eubacterial, human (including human mitochondrial) and fungal rRNA molecules.

It is another aspect of the present invention to provide probes and probe sets which provide a basis for discriminating between Gram positive and Gram negative eubacteria.

Methods for detecting, identifying and enumerating bacteria in normally sterile body fluids vary with the type of sample and the suspected pathogen. No currently available method is optimal for the detection of all pathogens. Often a combination of methods must be used to increase the likelihood that the pathogen will be detected. All commonly used methods for detection of, for example, bacteremia or bacterial septicemia rely on the *in vitro* cultivation of microbes from clinical samples. Generally, a blood sample is drawn from a patient and incubated in a rich artificial culture medium and monitored for 1 to 14 days. During this time, the medium is examined or blindly sub-cultured (plated), or assayed chemically or isotopically for evidence of bacterial growth or fermentative processes. Clinicians generally draw two or three samples of 10 milliliters of blood which may yield as few as one to ten colony forming units of bacteria for a positive diagnosis. Following the isolation of individual colonies of bacteria on diagnostic solid media and/or by Gram-staining, presumptive identification of the bacteria (or fungus) is made.

All cultivation methods suffer a number of serious shortcomings, including the following:

- High material costs;
- Labor intensive;
- Technologists extensively handle dangerous bodily fluids;
- False positives due to handling;
- False negatives due to low viable cell numbers;
- False negatives due to fastidious media requirements of many potential pathogens; and
- Relatively long time to positive diagnosis and identification.

Because of the relatively long time required by current methods to achieve a diagnosis and because of the potentially life threatening nature of such infections, antimicrobial therapy often is begun empirically before the results of such tests can be known.

Therefore, it is another aspect of the present invention to provide nucleic acid probes which are broadly specific for all eubacteria and which preferably do not react with other eukaryotic pathogens, especially fungi, which may be present in sampled materials.

It is yet another aspect of the present invention to provide probes which may be used in a variety of assay systems which avoid many of the disadvantages associated with traditional, multi-day culturing techniques.

It is still another aspect of the present invention to provide probes that are capable of hybridizing to the ribosomal ribonucleic acid (rRNA) of the targeted eubacterial organisms under normal assay conditions.

While Kohne et al. (Biophysical Journal 8:1104-1118, 1968) discuss one method for preparing probes to rRNA sequences, they do not provide the teaching necessary to make broad-specificity eubacterial probes.

Pace and Campbell (Journal of Bacteriology 107:543-547, 1971) discuss the homology of ribosomal ribonucleic acids from diverse bacterial species and a hybridization method for quantitating such homology levels. They do not identify particular nucleic acid sequences shared by bacteria, but absent in eukaryotes. Woese (Microbiological Reviews 51:221-271, 1987) describes the breadth of the eubacteria, in terms of rRNA sequence, but does not indicate sequences of interest for complete bacterial inclusivity. These references, however, fail to relieve the deficiency of Kohne's teaching with respect to eubacterial probes and, in particular, do not provide eubacterial specific probes useful in assays for detecting eubacteria in clinical or other samples.

Giovannoni et al. (Journal of Bacteriology 170:720-726, 1988) describe a number of probes which are claimed to be useful for the identification of broad groups of eubacteria, archaeobacteria and eukaryotes. However, Giovannoni et al. do not disclose the probes of the present invention. Nor do they provide the teaching necessary to design such probes.

Hogan et al. (European patent publication WO 88/03957) describe a number of probes which are claimed to hybridize to a broad representation of eubacteria. However, Hogan et al. do not teach the probes of the present invention and also fail to relieve the deficiency of Kohne's teaching with respect to these probes.

Ribosomes are of profound importance to all organisms because they serve as the only means of translating genetic information into cellular proteins. A clear manifestation of this importance is the observation that all cells have ribosomes. Actively growing bacteria may have 20,000 or more ribosomes per cell. This makes ribosomes one of the most abundant macromolecular entities in a cell, and an attractive diagnostic assay target.

Ribosomes contain three distinct RNA molecules which in *Escherichia coli* are referred to as 5S, 16S and 23S rRNAs. These names historically are related to the size of the RNA molecules, as determined by their sedimentation rate. In actuality, however, ribosomal RNA molecules vary in size between organisms. Nonetheless, 5S, 16S, and 23S rRNA are commonly used as generic names for the homologous RNA molecules in any bacteria, and this convention will be continued herein. Discussion will be confined to 16S and 23S rRNAs.

As used herein, probe(s) refer to synthetic or biologically produced nucleic acids (DNA or RNA) which, by design or selection, contain specific nucleotide sequences that allow them to hybridize under defined predetermined stringencies, specifically (i.e., preferentially, see below - Hybridization) to target nucleic acid sequences. In addition to their hybridization properties, probes also may contain certain constituents that pertain to their proper or optimal functioning under particular assay conditions. For example, probes may be modified to improve their resistance to nuclease degradation (e.g. by end capping), to carry detection ligands (e.g. fluorescein, 32-Phosphorous, biotin, etc.), or to facilitate their capture onto a solid support (e.g., poly-deoxyadenosine "tails"). Such modifications are elaborations on the basic probe function which is its ability to usefully discriminate between target and non-target organisms in a hybridization assay.

Hybridization traditionally is understood as the process by which, under predetermined reaction conditions, two partially or completely complementary strands of nucleic acid are allowed to come together in an antiparallel fashion (one oriented 5' to 3', the other 3' to 5') to form a double-stranded nucleic acid with specific and stable hydrogen bonds. (Note that nucleic acids do have a polarity; that is, one end of a nucleic acid strand is chemically different from another. This is defined by the polarity of the chemical linkages through the asymmetric sugar moiety of the nucleotide components. The terms 5' and 3' specifically refer to the ribose sugar carbons which bear those names. Except in rare or unusual circumstances, nucleic acid strands do not associate through hydrogen bonding of the base moieties in a parallel fashion. This concept is well understood by those skilled in the art.)

The stringency of a particular set of hybridization conditions is defined by the base composition of the probe/target duplex, as well as by the level and geometry of mispairing between the two nucleic acids.

Stringency may also be governed by such reaction parameters as the concentration and type of ionic species present in the hybridization solution, the types and concentrations of denaturing agents present, and/or the temperature of hybridization. Generally, as hybridization conditions become more stringent, longer probes are preferred if stable hybrids are to be formed. As a corollary, the stringency of the conditions under which a hybridization is to take place (e. g., based on the type of assay to be performed) will dictate certain characteristics of the preferred probes to be employed. Such relationships are well understood and can be readily manipulated by those skilled in the art.

As a general matter, dependent upon probe length, such persons understand stringent conditions to mean approximately 35°C-65°C in a salt solution of approximately 0.9 molar.

Summary of the Invention

In accordance with the various principles and aspects of the present invention, there are provided nucleic acid probes and probe sets comprising deoxyribonucleic acid (DNA) or ribonucleic acid (RNA) sequences which hybridize, under specific conditions, to the ribosomal RNA molecules (rRNA), rRNA genes (rDNA), and certain amplification and *in vitro* transcription products thereof of eubacteria but which do not hybridize, under the same conditions, to the rRNA or rDNA of eukaryotic cells which may be present in test samples. In addition, certain of the probes and probe sets described herein may be used as primers for the specific amplification of eubacterial rRNA or rDNA sequences which may be present in a sample by such methods as the polymerase chain reaction (US 4,683,202) or transcriptional amplification systems (e.g. TAS, Kwok et al., 1989, Proceedings of the National Academy of Science 86:1173-1177).

The probes of the present invention are capable of hybridising to rDNA and rRNA of eubacteria, but not to rRNA or rDNA of mouse L cells, wheat germ, human blood or *Candida albicans*, and have a sequence that is complementary to, or homologous to, at least 90% of sequence comprising any 10 consecutive nucleotides within any one of 22 probes (which are designated by a four digit number and whose sequences are given later).

The probes of the present invention advantageously provide the basis for development of valuable nucleic acid hybridization assays for the specific detection of eubacteria in clinical samples such as blood, urine, cerebrospinal fluid, biopsy, synovial fluid, or other tissue or fluid samples from humans or animals. The probes also provide the basis for testing, for example in quality control, substances that are presumed sterile, e.g., pharmaceuticals. The probes described herein are specifically complimentary to certain highly conserved bacterial 23S or 16S rRNA sequences.

The detection of bacteria by nucleic acid hybridization constitutes enhanced performance capability compared to the available culture-dependent tests for several reasons including:

- a) increased sensitivity; i.e., the ability to detect said bacteria in a given sample more frequently;
- b) potentially significant reductions in assay cost due to the use of inexpensive reagents and reduced labor;
- c) accurate detection of even nutritionally fastidious strains of bacteria;
- d) faster results because such tests do not require the isolation of the target bacterium from the sample prior to testing;
- e) the ability to screen, in a batch mode, a large number of samples, and only culture those identified as "hybridization positive";
- f) potential detection of phagocytized organisms eliminating the need for multiple, punctuated blood samples in order to sample the cyclical "window" of viable organisms (which probably depends on host immunological cycles);
- g) some reduction of technologist handling of potentially infectious body fluids;
- h) the ability to detect very low numbers of targets by amplifying either the bacterial signal or target using *in vitro* nucleic acid amplification.

It has been discovered that other advantages incurred by directing the probes of the present invention against rRNA include the fact that the rRNAs detected constitute a significant component of cellular mass. Although estimates of cellular ribosome content vary, actively growing *Escherichia coli*, for example, may contain upwards of 50,000 ribosomes per cell, and therefore 50,000 copies of each of the rRNAs (present in a 1:1:1 stoichiometry in ribosomes). The abundance of ribosomes in other bacteria particularly under other, less favorable, metabolic conditions may be considerably lower. However, under any circumstances, rRNAs are among the most abundant cellular nucleic acids present in all cell types. In contrast, other potential cellular target molecules such as genes or RNA transcripts thereof, are less ideal since they are present in much lower abundance.

A further unexpected advantage is that the rRNAs (and the genes specifying them) appear not to be subject to lateral transfer between contemporary organisms. Thus, the rRNA primary structure provides an organism-specific molecular target, rather than a gene-specific target as would likely be the case, for example of a plasmid-borne gene or product thereof which may be subject to lateral transmission between contemporary organisms.

Additionally, the present invention provides probes to eubacterial rRNA target sequences which are sufficiently similar in most or all eubacteria tested that they can hybridize to the target region in such eubacteria. Advantageously, these same rRNA target sequences are sufficiently different in most non-eubacterial rRNAs that, under conditions where the probes hybridize to eubacterial rRNAs they do not hybridize to most non-eubacterial rRNAs. These probe characteristics are defined as inclusivity and exclusivity, respectively.

The discovery that probes could be generated with the extraordinary inclusivity and exclusivity characteristics of those of the present invention with respect to eubacteria was unpredictable and unexpected.

The various aspects of the invention for which protection is sought are featured in appending claims.

Brief Description of the Figures

Further understanding of the principles and aspects of the present invention may be aided by reference to the tables wherein:

Figure 1 - Shows an evolutionary "tree" of the major genetic "kingdoms" of life (Woese, 1987, Microbiological Reviews 51:221- 271). The branching patterns represent the mutational distances between the 16S rRNA sequences of the represented organism. Such comparisons readily distinguish the eubacteria from the archaeobacteria and eukaryotes.

Figure 2 - Shows a more detailed evolutionary tree of the eubacterial kingdom (ibid.). So far about 10 major divisions/phyla have been defined based on 16S rRNA sequence comparisons. Certain discriminations among eubacterial divisions can be important in a clinical context and certain of the probes of the present invention do exhibit preferential hybridization to one or more of the eubacterial divisions. Therefore, the test organisms listed in Tables 3, 4 and 5 are grouped according to the divisions shown in

Figure 2 so that significant patterns of hybridization may be most easily discerned.

Brief Description of the Tables

Table 1 - Shows alignment of the nucleotide sequences of the preferred 16S rRNA-targeted probes of the present invention with their target nucleotide sequences in *E. coli* 16S rRNA. Very extensive sequence comparison to some 350 aligned 16S and 18S rRNA sequences were performed during the development of the probes of the present invention. It simply is not practical to show this analysis in detail. However, a consensus sequence (CONS-90%) of highly conserved 16S rRNA nucleotide positions is provided as a summary of the patterns of nucleotide sequence variation discovered among representative eubacteria. A nucleotide on the CONS-90% line indicates that that nucleotide is found at the homologous position in 90% or greater of the eubacterial sequences inspected. Note that the probe target regions all correspond to clusters of high sequence conservation among the eubacterial 16S and 23S rRNA molecules.

Since the *E. coli* 16S and 23S rRNA sequences were among the first full rRNA sequences obtained, the assigned position numbers have become a convenient and commonly accepted standard for explicitly identifying the homologous regions in other rRNA sequences under consideration. In Table 1, the *E. coli* RNA (target) sequence is written 5' to 3'. Probe sequences are DNA and written 3' to 5', except for probes 1638, 1642 and 1643 which are designed to hybridize to the rRNA-complementary sequence rather than the rRNA itself. These latter probes have the same "sense" (i. e. polarity) as the rRNA and are written 5' to 3'.

Table 2 - Shows alignment of the nucleotide sequences of the preferred 23S rRNA-targeted probes of the present invention with their target nucleotide sequences in *E. coli* 23S rRNA. As in Table 1 the *E. coli* sequence numbering is used as a standard in order to identify the homologous probe target sequences in all 23S rRNAs. CONS-90% has the same meaning as in TABLE 1. For the 23S rRNA analyses only about 30 sequences were available. However, these represent most of the major eubacterial divisions shown in Figure 2. In the probe 1730 sequence, "R" = a 1:1 mixture of A and G at that position.

Table 3 - Exemplifies the inclusivity and exclusivity behavior of a number of the preferred 16S rRNA-targeted probes toward a representative sampling of eubacterial and non-eubacterial rRNAs in a dot blot hybridization assay.

Table 4 - Exemplifies the inclusivity and exclusivity behavior of a number of the preferred 23S rRNA-targeted probes toward a representative sampling of eubacterial and non-eubacterial rRNAs in a dot blot hybridization assay.

Table 5 - Exemplifies the inclusivity and exclusivity behavior of a number of additional preferred 16S and 23S rRNA-targeted probes toward a representative sampling of eubacterial and non-eubacterial rRNAs in a dot blot hybridization assay. These probes exhibit useful patterns of hybridization to specific subgroups of eubacteria - notably Gram positive and Gram negative bacteria.

Detailed Description of the Invention and Best Mode

Probe Development Strategy:

The first step taken in the development of the probes of the present invention involved identification of regions of 16S and 23S rRNA which potentially could serve as target sites for eubacteria specific nucleic acid probes. This entailed finding sites which are:

- 1) highly conserved (few nucleotide changes, deletions, or insertions) among eubacterial rRNA sequences, and
- 2) substantially different in non-eubacterial rRNA sequences.

For this analysis, precise alignments of available 16S and 23S rRNA sequences were developed. A number of 16S and 23S rRNA sequences were determined as part of this effort. Such nucleotide sequences were determined by standard laboratory protocols either by cloning (Maniatis et al., 1982, Molecular Cloning; A Laboratory Manual, Cold Spring Harbor Laboratory, New York, pp 545) and sequencing (Maxam and Gilbert, 1977, Proceedings of the National Academy of Science, USA 74:560-564; Sanger et al., 1977, Proceedings of the National Academy of Science, USA 74:5463-5467) the genes which specify the rRNAs, and/or by direct sequencing of the rRNAs themselves using reverse transcriptase (Lane et al., 1985, Proceedings of the National Academy of Science, USA 82:6955-6959; Lane, manuscript in preparation).

A computer algorithm, operating on the aligned set of 16S and 23S rRNA sequences, was used to identify regions of greatest similarity among eubacteria. Nucleic acid probes to such regions will hybridize most widely among diverse eubacteria.

Such regions of homology among eubacteria next were assessed for differences with non-eubacterial rRNA sequences. In particular, sequence differences between eubacterial and human, fungal, and mitochondrial sequences were sought.

Forty one probes were designed based on these analyses; 22 targeting 23S rRNA and 19 targeting 16S rRNA.

The hybridization behavior of these probes toward extensive panels of eubacteria was determined by hybridization analysis in a dot blot format.

Physical Description of the Probes:

The foregoing probe selection strategy yielded a number of probes useful for identifying eubacteria in samples and include the following preferred oligonucleotide probes:

Probes within the invention are set out in the appended claims, and are nos: 1661, 1739, 1746, 1743, 1639, 1640, 1641, 1656, 1657, 1653, 1654, 1655, 1651, 1652, 1595, 1600, 1601, 1602, 1598, 1599, 1596 or 1597. All other probes are provided for comparison.

16S rRNA-targeted probes:

Probe 1638: 5'-AGAGTTTGATCCTGGCTCAG-3'

Probe 1642: 5'-AGAGTTTGATCATGGCTCAG-3'

Probe 1643: 5'-AGAGTTTGATCCTGGCTTAG-3'

Probe 1738: 5'-CTGAGCCAGGATCAAACCTCT-3'

Probe 1744: 5'-CAGCGTTCGTCCTGAGCCAGGATCAAACCT-3'

Probe 1659: 5'-CTGCTGCCTCCCGTAGGAGT-3'

Probe 1660: 5'-CTGCTGCCTCCCGTAGGAGTTTGGGCGGTGTCTCAGTTCCAGTGT-3'

*Probe 1661: 5'-TATTACCGCGGCTGCTGGCACGGAGTTAGCCG-3'

*Probe 1739: 5'-GCGTGGACTACCGGGGTATCTAATCCTGTTTGCTCCCCACGCTTTCG-3'

Probe 1740: 5'-GGGTTGCGCTCGTTGCGGGACTTAACCCGACATCTCACGGCACCAGCT
CACGACAGCCATGCAT-3'

- 5 Probe 1741: 5'-CTCACGGCAGGAGCTGACGACAGCCATGCAT-3'
Probe 1742: 5'-GGGTTGCGCTCGTTGCGGGACTTAACCCGACAT-3'
Probe 1745: 5'-AGCTGACGACAACCATGCACCACCTGT-3'
*Probe 1746: 5'-TCATAAGGGGCGATGATGATTTGACGTCAT-3'
*Probe 1743: 5'-GTACAAGGCCCGGGAACGTATTCACCG-3'
- 10 Probe 1637: 5'-AAGGAGGTGATCCAGCC-3'
*Probe 1639: 5'-ACGGTTACCTTGTTACGACTT-3'
*Probe 1640: 5'-ACGGCTACCTTGTTACGACTT-3'
*Probe 1641: 5'-ACGGATACCTTGTTACGACTT-3'
- 15 23S rRNA-targeted probes:
- Probe 1730: 5'-CTTTTCTCCTTTCCCTCRCGGTAAGTTCRGTATCGGTC-3'
Probe 1731: 5'-CTTTTTCGCCTTTCCCTCGCGGTAAGTTCGCTATCGGTC-3'
Probe 1658: 5'-TCTTTAAAGGGTGGCTGCTTCTAAGCCAACATCCTGGTTG-3'
- 20 *Probe 1656: 5'-CTACCTGTGTGCGTTTTCGCGTACGGGC-3'
*Probe 1657: 5'-GGTATTCTCTACCTGACCACCTGTGTCGGTTTGGGGTACG-3'
*Probe 1653: 5'-CCTTCTCCCGAAGTTACGGGGGCATTTTGCCTAGTTCCTT-3'
*Probe 1654: 5'-CCTTCTCCCGAAGTTACGGGGTTCATTTTGCCGAGTTCCTT-3'
*Probe 1655: 5'-CCTTCTCCCGAAGTTACGGCACCATTTTGCCGAGTTCCTT-3'
- 25 *Probe 1651: 5'-CTCCTCTTAACCTTCCAGCACCGGGCAGGC-3'
*Probe 1652: 5'-TTCGATCAGGGGCTTCGCTTTCGCTGACCCCATCAATTAA-3'
Probe 1512: 5'-TTAGGACCGTTATAGTTACGGCCGCGTTTACTGGGGCTT-3'
Probe 1256: 5'-GGTCGGAACCTTACCCGACAAGGAATTTGCTACCTTAG-3'
Probe 1398: 5'-GGTCGGTATTTAACCAGACAAGGAATTTGCTACCTTAG-3'
- 30 Probe 1511: 5'-CGTGCGGGTCGGAACCTTACCCGACAAGGAATTTGCTACCTAG-3'
*Probe 1595: 5'-CGATATGAACCTTTGGGCGGTATCAGCCTGTTATCCCCGG-3'
*Probe 1600: 5'-CAGCCCCAGGATGAGATGAGCCGACATCGAGGTGCCAAAC-3'
*Probe 1601: 5'-CAGCCCCAGGATGTGATGAGCCGACATCGAGGTGCCAAAC-3'
*Probe 1602: 5'-CAGCCCCAGGATGCGATGAGCCGACATCGAGGTGCCAAAC-3'
- 35 *Probe 1598: 5'-CGTACCGCTTTAAATGGCGAACAGCCATACCCTTGGGACC-3'
*Probe 1599: 5'-CGTGCCGCTTTAATGGCGAACAGCCCAACCCTTGGGACC-3'
*Probe 1596: 5'-GATAGGGACCGAAGTGTCTCACGACGTTTGAACCCAGCT-3'
*Probe 1597: 5'-GATAGGGACCGAAGTGTCTCACGACGTTCTGAACCCAGCT-3'

The specific behaviors of the aforementioned probes are dependent to a significant extent on the assay format in which they are employed. Conversely, the assay format will dictate certain of the optimal design features of particular probes. The "essence" of the probes of the invention is not to be construed as restricted to the specific string of nucleotides in the named probes. For example, the length of these particular oligonucleotides was optimized for use in the dot blot assay (and certain other anticipated assays) described below. It is well known to those skilled in the art that optimal probe length will be a function of the stringency of the hybridization conditions chosen and hence the length of the instant probes may be altered in accordance therewith. Also, in considering sets comprised of more than one probe, it is desirable that all probes behave in a compatible manner in any particular format in which they are employed. Thus, the exact length of a particular probe will to a certain extent, reflect its specific intended use. Again, given the probes of the instant invention, these are familiar considerations to one of ordinary skill in the art.

50 The "essence" of the probes described herein resides in the discovery and utilization of the specific sequences described above and given in Table 1 and Table 2.

* referred to in claim 1.

55

* referred to in claim 1.

* referred to in claim 1.

Hybridization Analysis of Probe Behavior:

The sequence comparisons which led to the discovery of the disclosed target sequences suggested that many of the probes should hybridize to a significant number of eubacteria. For the 16S rRNA analyses, some 350 sequences were considered in designing the probes; for the 23S rRNA analyses only about 30 eubacterial sequences were available. Since it is impossible to test every eubacterial strain, greater sequence variation might exist in other eubacterial strains not inspected by sequence analysis which might reduce or eliminate hybridization by the prospective probes to such untested eubacteria. As can be seen in Tables 3, 4 and 5, some probes of extremely broad inclusivity nevertheless fail to hybridize to certain bacteria. Therefore, carefully documenting the hybridization behavior to a large and representative number of eubacteria is an important element in documenting that such probes are capable of detecting all eubacteria or, failing that, for documenting which eubacteria are not detected. Such "failures" may not be clinically significant or alternatively, may be compensated for by appropriate inclusion of other probes of the instant inventions.

Equally as important as the inclusivity behavior of the probes, is their exclusivity behavior, i.e., their reactivity toward non-eubacteria. As mentioned, demonstrating a lack of hybridization to human and fungal rRNAs is of paramount importance in the types of clinical applications envisioned for such probes. Therefore, the behavior of the probes toward representative eubacterial, human and fungal rRNAs was determined by hybridization analysis using a dot blot procedure.

Example 1: Dot-blot analysis of probe hybridization behavior.

Dot-blot analysis, in accordance with well known procedures, involves immobilizing a nucleic acid or a population of nucleic acids on a filter such as nitrocellulose, nylon, or other derivatized membranes which readily can be obtained commercially, specifically for this purpose. Either DNA or RNA can be easily immobilized on such a filter and subsequently can be probed or tested for hybridization under any of a variety of conditions (i.e., stringencies) with nucleotide sequences or probes of interest. Under stringent conditions, probes whose nucleotide sequences have greater complementarity to the target sequence will exhibit a higher level of hybridization than probes containing less complementarity. For most of the oligonucleotide probes described herein, hybridization to rRNA targets at 60°C for 14-16 hours (in a hybridization solution containing 0.9 M NaCl, 0.12 M Tris-HCl, pH 7.8, 6 mM EDTA, 0.1 M KPO₄, 0.1% SDS, 0.1% pyrophosphate, 0.002% ficoll, 0.02% BSA, and 0.002% polyvinylpyrrolidone), followed by standard post-hybridization washes to remove unbound and non-specifically hybridized probe (at 60°C in 0.03 M NaCl, 0.004 M Tris-HCl, pH 7.8, 0.2 mM EDTA, and 0.1% SDS), would be sufficiently stringent to produce the levels of specificity demonstrated in Tables 3, 4 and 5. The exceptions to these conditions are probe 1738 (which was hybridized at 37°C), and probe 1746 (which was hybridized at 37°C and washed at 50°C).

Techniques also are available in which DNA or RNA present in crude (unpurified) cell lysates can be immobilized without first having to purify the nucleic acid in question (e.g. Maniatis, T., Fritsch, E. F. and Sambrook, J., 1982, Molecular Cloning: A Laboratory Manual).

The dot-blot hybridization data shown in Tables 3, 4 and 5 were generated by hybridization of the indicated probes to purified RNA preparations from the indicated eubacterial, fungal and human specimens. Bacterial and fungal RNAs were purified from pure cultures of the indicated organisms. Mouse RNA was purified from L cells (a tissue culture cell line). Wheat germ RNA was purified from a commercial preparation of that cereal product. Human blood and stool RNAs were purified from appropriate specimens obtained from normal, healthy individuals.

Purified RNA was used, rather than cell lysates for a number of simple technical reasons. The most important of these relate to proper interpretation of the relative signal arising from the hybridization of any particular probe to individual organisms. RNA content per cell is known to vary widely among different bacteria and varies even more between bacteria and eukaryotic cells. In addition, the specific metabolic status of cells at the time of harvest can have a profound influence on the amount and integrity of the RNA recovered. Some bacteria, for example, begin to degrade their RNA very rapidly upon reaching the stationary growth phase. The organisms represented in Tables 3, 4 and 5 comprise an extremely diverse collection in every respect. Represented are Gram positive and Gram negative bacteria, photosynthetic and chemosynthetic, heterotrophic and lithotrophic, and anaerobic and aerobic metabolisms. By using known, equivalent amounts of purified RNA in the individual "dots," relative levels of hybridization of each probe to each organism can be meaningfully compared without regard to the idiosyncracies of nucleic acid preparation from individual types of bacteria represented.

RNA was prepared by a variation on standard published methods which has been developed in our laboratory (W. Weisburg, unpublished). The method rapidly yields bulk high molecular weight RNA in a highly purified but relatively unfractionated form. Little or no DNA, or low molecular weight RNA species are found in RNA prepared in this fashion.

A large proportion of the RNA is 16S and 23S rRNA (18S and 28S rRNA in eukaryotes) as is true of the RNA in the intact cells. The method is rapid and convenient, but otherwise is not relevant to interpretation of the dot-blot results presented in Tables 3, 4 and 5. Most other currently accepted methods available in the literature which yield RNA of reasonable intactness will yield equivalent hybridization results.

For the hybridization experiments reported in Tables 3, 4 and 5, probes were end-labeled with radioactive 32-phosphorous, using standard procedures. Following hybridization and washing as described above, the hybridization filters were exposed to X-ray film and the intensity of the signal evaluated with respect to that of control RNA spots containing known amount of target RNA of known sequence.

A scale of hybridization intensity ranging from + + + + (hybridization signal equivalent to that of control spots) to + (barely detectable even after long exposure of the x-ray film) has been used to compare hybridization signals between organisms and probes. + + + signal indicates a very strong signal only slightly less intense than control spots. + + indicates a clearly discernible hybridization signal, but one that is noticeably weaker than the control spots. Note that while more "quantitative" ways to record hybridization signal are available, they are much more cumbersome to employ and, in our experience, not really any more useful for probe evaluation than the method employed in Tables 3, 4 and 5. In fact, because of certain uncontrollable variables in spotting exactly equivalent amounts of target RNA (of equivalent intactness) from such disparate organisms, numerically more precise counting methods are only deceptively more quantitative. In our experience, an organism generating a + + or greater signal to a particular probe is easily distinguished from one generating a "-" signal. This is true of a variety of assay formats that have been tested.


As is evident in Tables 3 and 4, 23S rRNA-targeted probes 1600, 1602, 1596, 1256 and 1512 and 16S rRNA-targeted probes, 1738, 1660, 1639, 1739, 1740, 1741 and 1743 hybridize most extensively among the eubacteria and are thus the most preferred. Other probes hybridize in a variety of patterns to subgroups of eubacteria and would be preferred for the detection of those subgroups or as components of more broadly inclusive probe sets. For example, probes 1599, 1656, 1744, 1745 and 1746 hybridize preferentially to Gram positive bacteria. Probes 1657, 1598 and 1595 hybridize preferentially to gram-negative bacteria, particularly to members of the so-called "purple bacterial" division (Figure 2 and Table 5).

Other probes exhibit other useful patterns of hybridization as is evident upon inspection of the data in Tables 3, 4 and 5. These probes can be combined in a variety of ways to create probe sets which exhibit the combined hybridization properties of the component probes. An example of one such hybridization format is given below (Example 2).

Alternatively, the probes could be used in a variety of subtractive hybridization schemes in which specific rRNA molecules are removed from the pool present in a mixed population of organisms prior to or simultaneous with the target organism-specific probes (e.g. Collins, European Patent Application 87309308.2, EP-A-0265244).

Example 2: Dual Probe, Sandwich Hybridization Assay

The probes of the present invention or derivatives thereof can be advantageously employed in a variety of other hybridization formats. One such format is a dual probe, sandwich-type hybridization assay such as that described, for example, in USSN 277,579; USSN 169,646, or USSN 233,683. In such a dual probe application, one probe (for example, probe 1602 or a derivative) would be ideally modified at its 3' terminus to contain a tract of about 20 - 200 deoxyadenosine (dA) residues. This would be used to "capture" the target rRNA (following liquid hybridization) from the test sample onto a solid support (e.g., beads, plastic surface, filter, etc.) which had been suitably derivatized with poly-deoxythymidine (dT) for this purpose. A second probe (for example, probe 1596 or derivative) would then be advantageously used as the detection probe and would be suitably derivatized with some detectable ligand (e.g. 32-P, fluorescein, biotin, etc.). Detection of the presence of the target nucleic acid in a test sample then would be indicated by capture of the detection ligand onto the solid surface through the series of hybridization interactions:

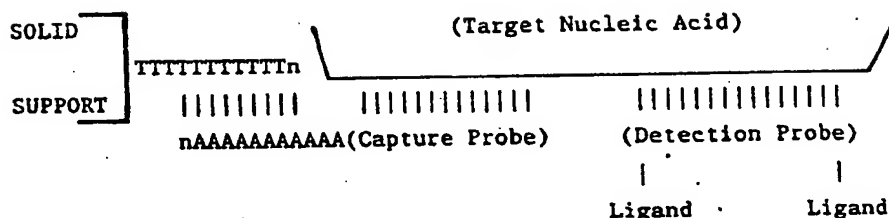

This could occur only if the target nucleic acid is present in the test sample. In principle, the above scheme could be employed with multiple capture and detection probes (probe sets) for the purpose of, for example, improving inclusivity or enhancing sensitivity of the assay.

Example 3: PCR Amplification of 16S rRNAs.

The polymerase chain reaction (PCR) is a well known method for amplifying target nucleic acid by "copying" the nucleic acid sequences located between two target sequences (US 4,683,202). The PCR process could be useful in an assay for the diagnosis of, for example, a non-viral pathogen by amplifying the genes encoding the pathogen's rRNA or rRNA genes and subsequently detecting that product. Implementation of this diagnostic strategy requires the invention of primers capable of amplifying the rRNA of the targeted organism(s). A second important application of such primers is in cloning amplified rRNA genes, and a third application is the direct sequencing of amplified rRNA genes.

Probes 1638, 1642, 1643, 1637, 1639*, 1640* and 1641* may be ideally used as primers for enzymatically copying and/or amplifying eubacterial 16S rRNAs or the genes encoding them. Details of the PCR procedure vary slightly depending on whether the target nucleic acid is single or double stranded, and whether it is DNA or RNA. However, the principle is the same in either case. Briefly, the steps are as follows:

- 1) Double-stranded DNA is denatured,
- 2) Oligonucleotide primers complementary to each of the sister DNA strands are annealed, and
- 3) deoxynucleotide triphosphate precursors are incorporated into newly synthesized sister DNA strands by extension of the primers from their 3' termini using DNA polymerase and/or reverse transcriptase.

Thus, a pair of oligonucleotide primers are required for the PCR reaction, one complementary to each strand within the target gene. They are positioned such that the newly synthesized product of one primer is also a target/template for the other primer. Thus the target nucleotide sequence located between the two primer annealing sites may be amplified many fold by repeating the steps listed above 20 to 30 times.

Probes 1638, 1642, 1643, 1637, 1639*, 1640* and 1641* are suitable for use as primers for enzymatically copying and/or amplifying eubacterial 16S rRNAs or the genes encoding them. That is, as a set, they will anneal very broadly among eubacterial rRNAs and rRNA genes and so will amplify any eubacterial rRNA sequences present in a sample.

Probes 1637, 1639*, 1640* and 1641* hybridize to the 16S rRNA (or rRNA-like strand of the ribosomal RNA gene) near its 3' end (Table 1). The template strand is read in the 3' to 5' direction producing an rRNA-complementary strand with the primer itself incorporated at its 5' terminus.

Probes 1638, 1642, and 1643 hybridize near the 5' end of the rRNA-complementary strand of the rRNA gene or to such a complement produced as described immediately above.

Individually, the above-described 16S rRNA amplification primers have approximately the following specificities:

5' primers:

- Probe 1638: most eubacteria
- Probe 1642: enterics and relatives
- Probe 1643: *Borrelia* spirochetes

* referred to in claim 1.

3' primers:

Probe 1637: most eubacteria

*Probe 1639: enterics, Deinococcus, Campylobacter

5 *Probe 1640: most eubacteria

*Probe 1641: fusobacteria, some Bacillus species

In test samples where the target bacterium is known, specific primers can be used. Where the target organism is not specifically known (for example, any eubacterium) all of the above mentioned primers can be used as a set.

10 The above described primers have been designed to amplify nearly the entire 16S rRNA sequence. Any of the other probes of the present invention or derivatives thereof can be used to amplify sub-segments of the 16S and 23S rRNAs or genes in a fashion similar to that just described.

Any such primers can be modified in a great number of ways to, for example, incorporate RNA polymerase promoters, cloning sites, etc. into the amplified transcripts.

15 While the description of the invention has been made with reference to detecting rRNA, it will be readily understood that the probes described herein and probes complementary to those described herein also will be useful for the detection of the genes (DNA) which specify the rRNA and, accordingly, such probes are encompassed within the present invention.

20

25

30

35

40

45

50

55

TABLE 1: 16S rRNA-TARGETED PROBES AND TARGET SEQUENCES

E. coli#s	8	38	
CONS-90%			
E. coli	AGAUUUGAUC UGGUCAG	GAACGCGGCG	
Probe 1738	5'-AAUUGAAGUUUGAUCUUGGUCAGAUUGAACGCGGCG		
Probe 1744	3'-TCCAAACTAGGACCGAGTC-5'		
Probe 1638	3'-TCAAACTAGGACCGAGTCCTGCTTGGAC-5'		
Probe 1642	5'-ACAGTTTGATCCTGGCTCAG-3'		
Probe 1643	5'-ACAGTTTGATCCTGGCTTAG-3'		
E. coli#s	313	357	
CONS-90%			
E. coli	G CACA UGG ACUGAGACACGG CCA ACUCCUACGGGAGGCGACAGU		
Probe 1660	GCCACACUGGAACUGAGACACGGUCCAGACUCCUACGGGAGGCGACAGUGG		
Probe 1659	3'-TGTGACCTTGACTCTGTGCGGGTTGAGGATGCCCTCCCTCTC-5'		
		3'-TCAGGATGCCCTCCCTCTC-5'	
E. coli#s	504	535	
CONS-90%			
E. coli	C CGGCUAACU GUGCCAGCAGCGCGGCUAAUAC		
Probe 1661	CACCGCUAACUCCGUGCCAGCAGCGCGGCUAAUACGG		
	3'-GCCGATTGAGGCACGGTCGTGCGCGGCATTAT-5'		
E. coli#s	764	810	
CONS-90%			
E. coli	G CGAAGCGUGGGAGC AACAGAUAGAUACCCUGGUAGUCCAGCG U		
Probe 1739	GUGCGAAGCGUGGGAGCAACAGAUAGAUACCCUGGUAGUCCAGCGCGU		
	3'-GCTTTCGCCACCCCTCGTTGCTTCTATCTATGGGGCCATCAGGTGCG-5'		

* referred to in claim 1.

TABLE 1 (cont'd): 16S rRNA-TARGETED PROBES AND TARGET SEQUENCES

E. coli#s	1044	1114
CONS-90%		
E. coli	ACAGGUG UGCAUGG UGUGUCAGGUCGUG CGUGAG UGUUGGGUUAAGUCCCGCAACGCCAACCC	
Probe 1745	GAGACAGGUGGUGCAUGGUGGUGUCAGGUCGUGUGUGAUAUGUUGGGUUAAGUCCCGCAACGCCAACCCUUA	
Probe 1740	3'-TGTCACCACTTACCAACAGCAGCTCGA-5'	
Probe 1741	3'-TAGCTACCGACAGCAGCTCGAGCAGCGGCGGTGCTCGCGGTGGG-5'	
Probe 1742	3'-TAGGTACCGACAGCAGCTCGAGCAGCGGCGGTGCTCGCGGTGGG-5'	
E. coli#s	1188	1216
CONS-90%		
E. coli	GGA GACGUCAA UC UCAUG CCCUUA G	
Probe 1746	GGGAUGACGUCAAUGCAUCAUGGCCCUUACGACCA	
	3'-TACTGCATTAGTAGTAGTACGGGGGAATACT-5'	
E. coli#s	1369	1395
CONS-90%		
E. coli	CGGUGAAUAGGUUC CGGG CUUGUACACA	
Probe 1743	CCACGGUGAAUAGGUUCCCGGGCCUUGUACACA	
	3'-GCCACTTATGCAAGGGCCCCCGAACATG-5'	
E. coli#s	1492	1541
CONS-90%		
E. coli	AAGUCGUAAACAAGGUA CC UA GAA UG GG UGGAU ACCUCCUUU	
Probe 1637	GUGAAGUCGUAAACAAGGUACCGUAGGGGAACCUCCGUUGGAUCACCUCCUUA-3'	
Probe 1639	3'-TTCAGCATGTTCATTCGCA-5'	
Probe 1640	3'-TTCAGCATGTTCATTCGCA-5'	
Probe 1641	3'-TTCAGCATGTTCATTCGCA-5'	

TABLE 2: 23S rRNA-TARGETED PROBES AND TARGET SEQUENCES

5	E. coli #s	442	481
	CONS-90%	GACCGAUAG G AC AGUACCGUGAGGGAAAGG GAAAAG AC	
	E.coli 23S	ACUGACCGAUAGUGAACCCAGUACCGUGAGGGAAAGCGAAAAGAAC	
	Probe 1730	3'-CTGGCTATCRCTTGGTCATGGCRCTCCCTTTCCRCTTTTC-5'	
10	Probe 1731	3'-CTGGCTATCGCTTGGTCATGGCGCTCCCTTTCCGCTTTTC-5'	
	E. coli #s	1049	1088
	CONS-90%	A ACA C AGGA GUUGGCUUAGAAGCAGCCA C UU AAAGA G	
15	E.coli 23S	AGACAGCCAGGAUGUUGGCUUAGAAGCAGCCAUCAUUUAAGAAAG	
	Probe 1658	3'-GTTGGTCCTACAACCGAATCTTCGTCCGTGGGAAATTTCT-5'	
	E. coli #s	1597	1639
	CONS-90%	CGUACC AAACCGACACAGGU G G A A C AG	
20	E.coli 23S	UCAAUUCGUACCCCAAACCGACACAGGUGGUCAGGUAGAGAAUACCAAG	
	Probe 1656	3'-CGGGCATGGCGTTTGGCTGTGTCCATC-5'	
	Probe 1657	3'-GCATGGGGTTTGGCTGTGTCCACCACTCATCTCTTATGG-5'	
25	E. coli #s	1664	1703
	CONS-90%	AAGGAACU GCAA U CCGUACUUCGG A AAGG	
	E.coli 23S	GUGAAGGAACUAGGC AAAAUGGUGCCGUAACUUCGGGAGAAGGCAC	
30	Probe 1653	3'-TTCCTTGATCCGTTTTACGGGGGCATTGAAGCCCTCTTCC-5'	
	Probe 1654	3'-TTCCTTGAGCCGTTTTACTGGGGCATTGAAGCCCTCTTCC-5'	
	Probe 1655	3'-TTCCTTGAGCCGTTTTACCACGGCATTGAAGCCCTCTTCC-5'	
	E. coli #s	1831	1860
35	CONS-90%	GAC CCUGCCC GUGC GGAAGGUUAA G	
	E.coli 23S	GACGCCUGCCCGGUGCCGGAAGGUUAAUUGAUGGGG	
	Probe 1651	3'-CGGACGGGCCACGACCTTCCAATTCCTC-5'	
40	E. coli #s	1851	1890
	CONS-90%	AGGUUAA G U G AAG A GAAGCC	
	E.coli 23S	AGGUUAAUUGAUGCCGUUAGCGCAAGCGAAGCUCUUGAUGGAAGCC	
45	Probe 1652	3'-AATTAAC TACCCCACTCGCGTTTCGCTTCGGGACTAGCTT-5'	
50			
55			

TABLE 2 (cont'd): 23S rRNA-TARGETED PROBES
AND TARGET SEQUENCES

5	E. coli #s	1889	1928
	CONS-90%	GAAGCCCC GU AACGGCGGCGCGUAACUAUAACGGUCCUAAGGU	
	E.coli 23S	UCGAAGCCCCCGGUAACGGCGGCGCGUAACUAUAACGGUCCUAAGGU	
10	Probe 1512	3'-TTCGGGGTCATTGCGCGCGGCATTGATATTGCCAGGATT-5'	
	E. coli #s	1925	1968
	CONS-90%	GUCCUAAGGUAGCGAAAUUCCUUGUCGGGUAAGUCCGACC GCACGAA	
15	E.coli 23S	GUCCUAAGGUAGCGAAAUUCCUUGUCGGGUAAGUCCGACCUGCACGAAU	
	Probe 1256	3'-GATTCCATCGCTTTAAGGAACAGCCCATTCAAGGCTGG-5'	
	Probe 1511	3'-CCATCGCTTTAAGGAACAGCCCATTCAAGGCTGGGCGTGC-5'	
	Probe 1398	3'-GATTCCATCGCTTTAAGGAACAGCCAATTTATGGCTGG-5'	
20	E. coli #s	2442	2481
	CONS-90%	AC C GGGGAUAACAGGCU AU C CC AG GU CA AUCG CG	
	E.coli 23S	ACUCCGGGGAUAACAGGCUGAUACCGCCCAAGAGUUCAUAUCCGACG	
25	Probe 1595	3'-GGCCCCCTATTGTCGACTATGGCGGGTTCTCAAGTATAGC-5'	
	E. coli #s	2490	2529
	CONS-90%	GUUUGGCACCUCGAUGUCCGCUC UC CAUCCUGGGGCGUG AG	
30	E.coli 23S	GGUGUUUGGCACCUCGAUGUCCGCUCAUCAUCCUGGGGCGUGAAG	
	Probe 1600	3'-CAAACCGTGGAGCTACAGCCGAGTAGAGTAGGACCCCGAC-5'	
	Probe 1601	3'-CAAACCGTGGAGCTACAGCCGAGTAGTGTAGGACCCCGAC-5'	
	Probe 1602	3'-CAAACCGTGGAGCTACAGCCGAGTAGCGTAGGACCCCGAC-5'	
35	E. coli #s	2535	2574
	CONS-90%	G GUUCCCAAGGGU GGCUGUUCGCC UUAAG GG ACG GA	
	E.coli 23S	GUAGGUCCCAAGGGUAUGGCUGUUCGCCAUUUAAGUGGUACCGCA	
40	Probe 1598	3'-CCAGGGTTCCCATACCGACAAGCGGTAAATTTCCGCCATGC-5'	
	Probe 1599	3'-CCAGGGTTCCCAACCCGACAAGCGGGTAATTTCCGCCGTGC-5'	
	E. coli #s	2577	2616
	CONS-90%	G GAGCUGGGU A AACGUCGUGAGACAGU GGUC CUAUC	
45	E.coli 23S	CGGAGCUGGGUUAGAAGCUGGUGAGACAGUCCGUCCCUAUCUGC	
	Probe 1596	3'-TCGACCCAAGTTTTGCAGCACTCTGTCAAGCCAGGGATAG-5'	
	Probe 1597	3'-TCGACCCAAGTCTTGCAGCACTCTGTCAAGCCAGGGATAG-5'	
50			
55			

TABLE 3: DOT BLOT HYBRIDIZATION OF 16S rRNA-TARGETED PROBES

Genus species	strain	div	PROBE HYBRIDIZATION									
			1738	1739	1659	1660	1661	1740	1741	1742	1743	1743
Acinetobacter calcoaceticus	GT0002	Purple	+++	+++	+++	+++	+++	+++	+++	+++	+++	+++
Aeromonas sobria	GT0007	gamma	+++	+++	+++	+++	+++	+++	+++	+++	+++	+++
Altremonas putrefaciens	GT1945	"	+++	+++	+++	+++	+++	+++	+++	+++	+++	+++
Citrobacter ananonicus	GT0690	"	+++	+++	+++	+++	+++	+++	+++	+++	+++	+++
Citrobacter diversus	GT0030	"	+++	+++	+++	+++	+++	+++	+++	+++	+++	+++
Citrobacter freundii	GT0687	"	+++	+++	+++	+++	+++	+++	+++	+++	+++	+++
Edwardsiella tarda	GT0569	"	+++	+++	+++	+++	+++	+++	+++	+++	+++	+++
Enterobacter agglomerans	GT0683	"	+++	+++	+++	+++	+++	+++	+++	+++	+++	+++
Enterobacter cloacae	GT0686	"	+++	+++	+++	+++	+++	+++	+++	+++	+++	+++
Enterobacter sakazakii	GT0062	"	+++	+++	+++	+++	+++	+++	+++	+++	+++	+++
Escherichia coli	GT1665	"	+++	+++	+++	+++	+++	+++	+++	+++	+++	+++
Escherichia coli	GT1592	"	+++	+++	+++	+++	+++	+++	+++	+++	+++	+++
Escherichia coli	GT1659	"	+++	+++	+++	+++	+++	+++	+++	+++	+++	+++
Haemophilus influenza	GT0244	"	+++	+++	+++	+++	+++	+++	+++	+++	+++	+++
Haemophilus ducreyi	GT0243	"	+++	+++	+++	+++	+++	+++	+++	+++	+++	+++
Hafnia alvei	GT0241	"	+++	+++	+++	+++	+++	+++	+++	+++	+++	+++
Morganella morganii	GT0303	"	+++	+++	+++	+++	+++	+++	+++	+++	+++	+++
Klebsiella pneumoniae	GT1500	"	+++	+++	+++	+++	+++	+++	+++	+++	+++	+++
Proteus mirabilis	GT1496	"	+++	+++	+++	+++	+++	+++	+++	+++	+++	+++
Providencia alcalifaciens	GT0371	"	+++	+++	+++	+++	+++	+++	+++	+++	+++	+++
Pseudomonas aeruginosa	GT1909	"	+++	+++	+++	+++	+++	+++	+++	+++	+++	+++
Salmonella arizona	GT0799	"	+++	+++	+++	+++	+++	+++	+++	+++	+++	+++
Salmonella typhimurium	GT0389	"	+++	+++	+++	+++	+++	+++	+++	+++	+++	+++
Serratia marcescens	GT0392	"	+++	+++	+++	+++	+++	+++	+++	+++	+++	+++
Shigella flexneri	GT0798	"	+++	+++	+++	+++	+++	+++	+++	+++	+++	+++
Vibrio parahaemolyticus	GT0568	"	+++	+++	+++	+++	+++	+++	+++	+++	+++	+++
Xanthomonas maltophilia	GT0417	"	+++	+++	+++	+++	+++	+++	+++	+++	+++	+++
Yersinia enterocolitica	GT0419	"	+++	+++	+++	+++	+++	+++	+++	+++	+++	+++
Alcaligenes faecalis	GT0610	Purple	+++	+++	+++	+++	+++	+++	+++	+++	+++	+++
Branhamella catarrhalis	GT0014	beta	+++	+++	+++	+++	+++	+++	+++	+++	+++	+++
Chromobacterium violaceum	GT2022	"	+++	+++	+++	+++	+++	+++	+++	+++	+++	+++
Kingella indologenes	0246	"	+++	+++	+++	+++	+++	+++	+++	+++	+++	+++
Moraxella osloensis	GT0301	"	+++	+++	+++	+++	+++	+++	+++	+++	+++	+++
Moraxella phenylpyruvica	GT0302	"	+++	+++	+++	+++	+++	+++	+++	+++	+++	+++

TABLE 3 (cont'd): DOT BLOT HYBRIDIZATION OF 16S rRNA-TARGETED PROBES

Genus species	strain	div	PROBE HYBRIDIZATION									
			1738	1739	1659	1660	1661	1740	1741	1742	1743	
<i>Borrelia burgdorferi</i>		Spiro	+++	+++	+++	+++	+++	+++	+++	+++	+++	
<i>Borrelia turicatae</i>		"	+++	+++	+++	+++	+++	+++	+++	+++	+++	
<i>Leptospira interrogans-pomona</i>		"	+++	+	+++	+++	+++	+++	+++	+++	+	
<i>Leptospira biflexa (Patoc-Patoc)</i>		"	+++	+	+++	+++	+++	+++	+++	+++	+	
<i>Leptospira biflexa (CDC)</i>		"	+++	+++	+++	+++	+++	+++	+++	+++	+++	
<i>Spirochaeta aurantia</i>		"	+++	+++	+++	+++	+++	+++	+++	+++	+++	
<i>Bacteroides fragilis</i>	25285 Bact		+++	+	+++	+++	+++	+++	+++	+++	+++	
<i>Bacteroides fragilis</i>	29771 "		+++	+	+++	+++	+++	+++	+++	+++	+++	
<i>Bacteroides thetaiotaomicron</i>	0572 "		+++	+	+++	+++	+++	+++	+++	+++	+++	
<i>Bacteroides melaninogenicus</i>	0011 "		+++	+	+++	+++	+++	+++	+++	+++	+++	
<i>Flavobacterium meningosepticum</i>	0237 "		+++	+++	+++	+++	+++	+++	+++	+++	+++	
<i>Chlamydia psittaci</i>		Chlam	-	+++	-	+++	+++	+++	+++	+++	+++	
<i>Chlamydia trachomatis</i>		"	-	+++	-	+++	+++	+++	+++	+++	+++	
<i>Chlorobium limicola</i>	Y400	Misc	+++	+++	+++	+++	+++	+++	+++	+++	+	
<i>Chloroflexus aurantiacus</i>	2608	"	+++	+++	+++	+++	+++	+++	+++	+++	+++	
<i>Deinococcus radiodurans</i>	2577	"	+++	+++	-	+++	+++	+++	+++	+++	-	
<i>Planctomyces maris</i>			+++	+++	+++	+++	+++	+++	+++	+++	+++	
Normal Stool RNA			-	+	-	-	+	-	-	-	-	
Mouse L-Cell			-	-	-	-	+	-	-	-	-	
Wheat Germ			-	-	-	-	+	-	-	-	-	
Normal Human Blood			-	-	-	-	+	-	-	-	-	
<i>Candida lusitanae</i>	403-87		-	-	-	-	+	-	-	-	-	
<i>Candida parapsilosis</i>	882-88		-	-	-	-	+	-	-	-	-	
<i>Candida tropicalis</i>	224-87		-	-	-	-	+	-	-	-	-	
<i>Candida albicans</i>	1008-88		-	-	-	-	+	-	-	-	-	
<i>Candida albicans</i>	223-87		-	-	-	-	+	-	-	-	-	
<i>Candida albicans</i>	819-88		-	-	-	-	+	-	-	-	-	

* Inclusivity and Exclusivity data was determined after overnight exposures.

** Each organism is represented by 100ng of CsTFA purified RNA.

*** Probe 1738 - hybridizations and washes were carried out at 37 C. +++ = positive control level of hybridization, + = barely detectable and - = zero, ND = not done.

TABLE 3 (cont'd): DOT BLOT HYBRIDIZATION OF 16S rRNA-TARGETED PROBES

Genus species	strain	div	PROBE HYBRIDIZATION										
			1738	1739	1659	1660	1661	1740	1741	1742	1743		
Mycoplasma pneumoniae	ATCC15531	"	+++	+++	+++	+++	+++	+++	+++	+++	+++	+	+
Mycoplasma putrefaciens	ATCC15718	"	+++	+++	+++	+++	+++	+++	+++	+++	+++	+	+
Peptostreptococcus productus	ATCC27340	"	+++	+++	+++	+++	+++	+++	+++	+++	+++	+	+
Planococcus citreus	ATCC14404	"	+++	+++	+++	+++	+++	+++	+++	+++	+++	+	+
Staphylococcus aureus	GT0399	"	+++	+++	+++	+++	+++	+++	+++	+++	+++	+	+
Staphylococcus aureus	GT1711	"	+++	+++	+++	+++	+++	+++	+++	+++	+++	+	+
Staphylococcus epidermidis	GT0401	"	+++	+++	+++	+++	+++	+++	+++	+++	+++	+	+
Staphylococcus haemolyticus	ATCC29970	"	+++	+++	+++	+++	+++	+++	+++	+++	+++	+	+
Streptococcus agalactiae	GT0405	"	+++	+++	+++	+++	+++	+++	+++	+++	+++	+	+
Streptococcus bovis	GT0668	"	+++	+++	+++	+++	+++	+++	+++	+++	+++	+	+
Streptococcus faecalis	GT0406	"	+++	+++	+++	+++	+++	+++	+++	+++	+++	+	+
Streptococcus morbillorum	GT2194	"	++	+++	+++	+++	+++	+++	+++	+++	+++	+	+
Streptococcus mutans	GT0412	"	+++	+++	+++	+++	+++	+++	+++	+++	+++	+	+
Streptococcus pneumoniae	GT0408	"	+++	+++	+++	+++	+++	+++	+++	+++	+++	+	+
Streptococcus salivarius	GT0410	"	+++	+++	+++	+++	+++	+++	+++	+++	+++	+	+
Streptococcus sanguis	GT0411	"	+++	+++	+++	+++	+++	+++	+++	+++	+++	+	+
Bifidobacterium dentium	GT0012	hiG+C	-	+++	+++	+++	+++	+++	+++	+++	+++	+	+
Corynebacterium genitalium	GT0045	Gm +	+++	+++	+++	+++	+++	+++	+++	+++	+++	+	+
Corynebacterium glutamicum	GT2120	"	+++	+++	+++	+++	+++	+++	+++	+++	+++	+	+
Corynebacterium pseudodiphtheriticum	GT2119	"	+++	+++	+++	+++	+++	+++	+++	+++	+++	+	+
Corynebacterium pseudotuberculosis	GT2122	"	+++	+++	+++	+++	+++	+++	+++	+++	+++	+	+
Corynebacterium pyogenes	GT2121	"	++	+++	+++	+++	+++	+++	+++	+++	+++	+	+
Corynebacterium xerosis	GT0046	"	++	+++	+++	+++	+++	+++	+++	+++	+++	+	+
Corynebacterium bovis	RCG	"	++	+++	+++	+++	+++	+++	+++	+++	+++	+	+
Mycobacterium kansasii	GT2191	"	+++	+++	+++	+++	+++	+++	+++	+++	+++	+	+
Nocardia asteroides	Misc	"	+++	+++	+++	+++	+++	+++	+++	+++	+++	+	+
Rhodococcus rhodochrous	GT2116	Gm +	+++	+++	+++	+++	+++	+++	+++	+++	+++	+	+
Aerococcus viridans	GT0238	"	+++	+++	+++	+++	+++	+++	+++	+++	+++	+	+
Fusobacterium necrophorum	ATCC27768	"	+++	+++	+++	+++	+++	+++	+++	+++	+++	+	+
Fusobacterium prausnitzii	GT2118	"	+++	+++	+++	+++	+++	+++	+++	+++	+++	+	+
Geella haemolysans	"	"	+++	+++	+++	+++	+++	+++	+++	+++	+++	+	+
Helicobacillus mobilis	"	"	+++	+++	+++	+++	+++	+++	+++	+++	+++	+	+
Phoridium ectocarpus	Cyano	"	+++	++	+++	+++	+++	+++	+++	+++	+++	+	+
Plectonema boryanum	"	"	+++	++	+++	+++	+++	+++	+++	+++	+++	+	+

TABLE 3 (cont'd): DOT BLOT HYBRIDIZATION OF 16S rRNA-TARGETED PROBES

			PROBE HYBRIDIZATION									
Genus species	strain	div	1738	1739	1659	1660	1661	1740	1741	1742	1743	
<i>Neisseria gonorrhoeae</i>	GT0315	.	+++	+++	+++	+++	+++	+++	+++	+++	+++	
<i>Neisseria meningitidis</i>	GT0349	.	+++	+++	+++	+++	+++	+++	+++	+++	+++	
<i>Pseudomonas acidovorans</i>	GT0376	.	+++	+++	+++	+++	+++	+++	+++	+++	+++	
<i>Pseudomonas cepacia</i>	GT2015	.	+++	+++	+++	+++	+++	+++	+++	+++	+++	
<i>Rhodocyclus gelatinosa</i>	ATCC17013	.	+++	+++	+++	+++	+++	+++	+++	+++	+++	
<i>Vitreoscilla stercoraria</i>	GT0810	Purple	+++	+++	+++	+++	+++	+++	+++	+++	+++	
<i>Achromobacter xerosis</i>	ATCC33463	alpha	+++	+++	+++	+++	+++	+++	+++	+++	+++	
<i>Acidiphilium cryptum</i>	GT2021	.	+++	+++	+++	+++	+++	+++	+++	+++	+++	
<i>Agrobacterium tumefaciens</i>	ATCC23448	.	+++	+++	+++	+++	+++	+++	+++	+++	+++	
<i>Brucella abortus</i>	GT2025	.	+++	+++	+++	+++	+++	+++	+++	+++	+++	
<i>Flavobacterium capsulatum</i>	GT2023	.	+++	+++	+++	+++	+++	+++	+++	+++	+++	
<i>Mycoplana bullata</i>	GT2020	.	+++	+++	+++	+++	+++	+++	+++	+++	+++	
<i>Pseudomonas diminuta</i>	ATCC17023	.	+++	+++	+++	+++	+++	+++	+++	+++	+++	
<i>Rhodobacter sphaeroides</i>	ATCC25903	.	+++	+++	+++	+++	+++	+++	+++	+++	+++	
<i>Rhodospirillum rubrum</i>	ATCC25364	.	+++	+++	+++	+++	+++	+++	+++	+++	+++	
<i>Thiobacillus versutus</i>	ATCC 7757	delta	+++	+++	+++	+++	+++	+++	+++	+++	+++	
<i>Desulfovibrio desulfuricans</i>	GT2095	Gm -	+++	+++	+++	+++	+++	+++	+++	+++	+++	
<i>Cardiobacterium hominis</i>	GT2172	Gm -	+++	+++	+++	+++	+++	+++	+++	+++	+++	
<i>Francisella tularensis</i>	GT0022	Campy	+++	+++	+++	+++	+++	+++	+++	+++	+++	
<i>Campylobacter jejuni</i>	GT0026	.	+++	+++	+++	+++	+++	+++	+++	+++	+++	
<i>Campylobacter pylori</i>	GT0027	.	+++	+++	+++	+++	+++	+++	+++	+++	+++	
<i>Campylobacter sputorum</i>	GT0803	lowG+C	+++	+++	+++	+++	+++	+++	+++	+++	+++	
<i>Bacillus brevis</i>	GT0804	Gm +	+++	+++	+++	+++	+++	+++	+++	+++	+++	
<i>Bacillus subtilis</i>	ATCC25537	.	+++	+++	+++	+++	+++	+++	+++	+++	+++	
<i>Clostridium clostridioforme</i>	GT0567	.	+++	+++	+++	+++	+++	+++	+++	+++	+++	
<i>Clostridium sordellii</i>	ATCC 3587	.	+++	+++	+++	+++	+++	+++	+++	+++	+++	
<i>Clostridium sporogenes</i>	ATCC19401	.	+++	+++	+++	+++	+++	+++	+++	+++	+++	
<i>Clostridium histolyticum</i>	ATCC13124	.	+++	+++	+++	+++	+++	+++	+++	+++	+++	
<i>Clostridium perfringens</i>	ATCC25582	.	+++	+++	+++	+++	+++	+++	+++	+++	+++	
<i>Clostridium ramosum</i>	ATCC33403	.	+++	+++	+++	+++	+++	+++	+++	+++	+++	
<i>Kurthia zopfii</i>	GT0256	.	+++	+++	+++	+++	+++	+++	+++	+++	+++	
<i>Lactobacillus acidophilus</i>	GT0258	.	+++	+++	+++	+++	+++	+++	+++	+++	+++	
<i>Lactobacillus plantarum</i>	IG3299	.	+++	+++	+++	+++	+++	+++	+++	+++	+++	
<i>Listeria monocytogenes</i>			+++	+++	+++	+++	+++	+++	+++	+++	+++	

TABLE 4: DOT BLOT HYBRIDIZATION OF 23S rRNA-TARGETED PROBES

Genus species	strain	div	PROBE HYBRIDIZATION									
			1730	1731	1658	1653	1654	1655	1651			
Acinetobacter calcoaceticus	GT0002	Purple	+++	+++	+++	+++	+++	+++	+++			
Aeromonas sobria	GT0007	gamma	+++	+++	+++	+++	+++	+++	+++			
Alteromonas putrefaciens	GT1945	"	+++	+++	+++	+++	+++	+++	+++			
Citrobacter amalonaticus	GT0690	"	+++	+++	+++	+++	+++	+++	+++			
Citrobacter diversus	GT0030	"	+++	+++	+++	+++	+++	+++	+++			
Citrobacter freundii	GT0687	"	+++	+++	+++	+++	+++	+++	+++			
Edwardsiella tarda	GT0569	"	+++	+++	+++	+++	+++	+++	+++			
Enterobacter agglomerans	GT0683	"	+++	+++	+++	+++	+++	+++	+++			
Enterobacter cloacae	GT0686	"	+++	+++	+++	+++	+++	+++	+++			
Enterobacter sakazakii	GT0062	"	+++	+++	+++	+++	+++	+++	+++			
Escherichia coli	GT1665	"	+++	+++	+++	+++	+++	+++	+++			
Escherichia coli	GT1592	"	+++	+++	+++	+++	+++	+++	+++			
Escherichia coli	GT1659	"	+++	+++	+++	+++	+++	+++	+++			
Haemophilus influenza	GT0244	"	+++	+++	+++	+++	+++	+++	+++			
Haemophilus influenza	ATCC33391	"	ND	ND	ND	ND	ND	ND	ND			
Haemophilus ducreyi	GT0243	"	+++	+++	+++	+++	+++	+++	+++			
Hafnia alvei	GT0241	"	+++	+++	+++	+++	+++	+++	+++			
Morganella morganii	GT0303	"	+++	+++	+++	+++	+++	+++	+++			
Klebsiella pneumoniae	GT1500	"	+++	+++	+++	+++	+++	+++	+++			
Proteus mirabilis	GT1496	"	+++	+++	+++	+++	+++	+++	+++			
Providencia alcalifaciens	GT0371	"	+++	+++	+++	+++	+++	+++	+++			
Pseudomonas aeruginosa	GT1909	"	+++	+++	+++	+++	+++	+++	+++			
Salmonella arizona	GT0799	"	+++	+++	+++	+++	+++	+++	+++			
Salmonella typhimurium	GT0389	"	+++	+++	+++	+++	+++	+++	+++			
Serratia marcescens	GT0392	"	+++	+++	+++	+++	+++	+++	+++			
Shigella flexneri	GT0798	"	+++	+++	+++	+++	+++	+++	+++			
Vibrio parahaemolyticus	GT0568	"	+++	+++	+++	+++	+++	+++	+++			
Xanthomonas maltophilia	GT0417	"	+++	+++	+++	+++	+++	+++	+++			
Yersinia enterocolitica	GT0419	"	+++	+++	+++	+++	+++	+++	+++			
Alcaligenes faecalis	GT0610	Purple	+++	+++	+++	+++	+++	+++	+++			
Branhamella catarrhalis	GT0014	beta	+++	+++	+++	+++	+++	+++	+++			
Chromobacterium violaceum	GT2022	"	+++	+++	+++	+++	+++	+++	+++			
Kingella indologenes	0246	"	+++	+++	+++	+++	+++	+++	+++			
Moraxella osloensis	GT0301	"	+++	+++	+++	+++	+++	+++	+++			

TABLE 4 (cont'd): DOT BLOT HYBRIDIZATION OF 23S rRNA-TARGETED PROBES

Genus species	strain	div	PROBE HYBRIDIZATION							
			1730	1731	1658	1653	1654	1655	1651	
<i>Moraxella phenylpyruvica</i>	GT0302	"	+++	+++	+	+++	+++	+++	+++	
<i>Morexelia gonorrhoeae</i>	GT0315	"	+++	+++	+++	-	+++	+++	+	
<i>Neisseria gonorrhoeae</i>	GT0349	"	+++	+++	+++	-	+++	+++	+++	
<i>Neisseria meningitidis</i>	GT0376	"	+++	+++	+++	-	+++	+++	+++	
<i>Pseudomonas acidovorans</i>	GT2015	"	+++	+++	+++	-	+++	+++	-	
<i>Pseudomonas cepacia</i>	ATCC17013	"	+++	+++	+++	-	+++	+++	+++	
<i>Rhodocyclus gelatinosa</i>		"	+++	+++	+++	-	+++	+++	+++	
<i>Vitreoscilla stercoraria</i>	GT0810	Purple	+++	+++	+++	-	+++	+++	+++	
<i>Achromobacter xerosis</i>	ATCC33463	alpha	+	-	++	-	-	-	+++	
<i>Acidiphilium cryptum</i>	GT2021	"	+++	++	+++	-	-	-	+++	
<i>Agrobacterium tumefaciens</i>	ATCC23448	"	+++	+++	+++	-	-	-	+++	
<i>Brucella abortus</i>	GT2025	"	+++	+++	+++	-	-	-	+++	
<i>Flavobacterium capsulatum</i>	GT2023	"	+++	+++	+	-	-	-	-	
<i>Mycoplana bullata</i>	GT2020	"	+++	+++	+	-	-	-	-	
<i>Pseudomonas diminuta</i>	ATCC17023	"	+++	+++	+++	-	-	-	+++	
<i>Rhodobacter sphaeroides</i>	ATCC25903	"	+++	+++	+++	-	+++	+++	+++	
<i>Rhodospirillum rubrum</i>	ATCC25364	"	+++	+++	+++	-	-	-	+++	
<i>Thiobacillus versutus</i>	ATCC 7757	delta	+++	+	++	+	+++	+++	+	
<i>Desulfovibrio desulfuricans</i>	GT2095	Gm -	+++	+++	-	-	+++	+++	+	
<i>Cardiobacterium hominis</i>	GT2172	Gm -	+++	+++	+++	-	-	-	+	
<i>Francisella tularensis</i>	GT0022	Campy	+++	+	+++	-	-	-	-	
<i>Campylobacter jejuni</i>	GT0026	"	+++	+	+++	-	-	-	+++	
<i>Campylobacter pylori</i>	GT0027	"	+++	+	+++	-	-	-	-	
<i>Campylobacter sputorum</i>	GT0803	lowG+C	+++	+++	+++	+++	+++	+++	+++	
<i>Bacillus brevis</i>	GT0804	Gm +	+++	+++	+++	+++	+++	+++	+++	
<i>Bacillus subtilis</i>	ATCC25537	"	-	-	+	+++	+++	+++	+++	
<i>Clostridium clostridioforme</i>	GT0567	"	+++	+++	+++	+++	+++	+++	+++	
<i>Clostridium sordellii</i>	ATCC 3587	"	+++	+++	+++	+++	+++	+++	+++	
<i>Clostridium sporogenes</i>	ATCC19401	"	+++	+++	+++	+++	+++	+++	+++	
<i>Clostridium histolyticum</i>	ATCC13124	"	+++	+++	+++	-	+++	+++	+++	
<i>Clostridium perfringens</i>	ATCC25582	"	+++	+++	+++	-	-	-	+++	
<i>Clostridium ramosum</i>	ATCC33403	"	+++	+++	+++	+++	+++	+++	+++	
<i>Kurthia zopfii</i>	GT0256	"	+++	+++	-	+	+++	+++	+++	
<i>Lactobacillus acidophilus</i>	GT0258	"	+++	+++	+++	+++	+++	+++	+++	
<i>Lactobacillus plantarum</i>		"	+++	+++	+++	+++	+++	+++	+++	

TABLE 4 (cont'd): DOT BLOT HYBRIDIZATION OF 23S rRNA-TARGETED PROBES

Genus species	strain	div	PROBE HYBRIDIZATION									
			1730	1731	1658	1653	1654	1655	1651			
<i>Listeria monocytogenes</i>	IG3299	"	+++	+++	+++	+++	+++	+++	+++			
<i>Mycoplasma pneumoniae</i>	ATCC15531	"	+++	++	+++	-	-	-	-			
<i>Mycoplasma putrefaciens</i>	ATCC15718	"	+++	++	++	-	-	-	-			
<i>Peptostreptococcus productus</i>	ATCC27340	"	-	-	++	+++	+++	+++	+++			
<i>Planococcus citreus</i>	ATCC14404	"	+++	+++	+++	+++	+++	+++	+++			
<i>Staphylococcus aureus</i>	GT0399	"	+++	+++	+++	+++	+++	+++	+++			
<i>Staphylococcus aureus</i>	GT1711	"	+++	+++	+++	+++	+++	+++	+++			
<i>Staphylococcus epidermidis</i>	GT0401	"	+++	+++	+++	+++	+++	+++	+++			
<i>Staphylococcus haemolyticus</i>	ATCC29970	"	+++	+++	+++	+++	+++	+++	+++			
<i>Streptococcus agalactiae</i>	GT0405	"	+++	+++	-	+++	+++	+++	+++			
<i>Streptococcus bovis</i>	GT0668	"	+++	+++	-	+++	+++	+++	+++			
<i>Streptococcus faecalis</i>	GT0406	"	+++	+++	+++	+++	+++	+++	+++			
<i>Streptococcus morbillorum</i>	GT2194	"	+++	++	+++	+++	+++	+++	+++			
<i>Streptococcus mutans</i>	GT0412	"	+++	++	-	+++	+++	+++	+++			
<i>Streptococcus pneumoniae</i>	GT0408	"	+++	+++	-	+++	+++	+++	+++			
<i>Streptococcus salivarius</i>	GT0410	"	+++	++	-	+++	+++	+++	+++			
<i>Streptococcus sanguis</i>	GT0411	"	+++	+++	-	+++	+++	+++	+++			
<i>Bifidobacterium dentium</i>	GT0012	hiG+C	+++	+++	-	-	-	-	-			
<i>Corynebacterium genitalium</i>	GT0045	Gm +	+++	+++	+++	+++	+++	+++	+++			
<i>Corynebacterium glutamicum</i>	GT2120	"	+++	+++	+++	+++	+++	+++	+++			
<i>Corynebacterium pseudodiphtheriticum</i>	GT2119	"	+++	+++	+++	+++	+++	+++	+++			
<i>Corynebacterium pseudotuberculosis</i>	GT2122	"	+++	+++	-	+++	+++	+++	+++			
<i>Corynebacterium pyogenes</i>	GT2121	"	+	-	+	+++	+++	+++	+++			
<i>Corynebacterium xerosis</i>	GT0046	"	-	-	+	+++	+++	+++	+++			
<i>Mycobacterium bovis</i>	BCG	"	+	-	+	+++	+++	+++	+++			
<i>Mycobacterium kansasii</i>	GT2191	"	+++	+++	+++	+++	+++	+++	+++			
<i>Nocardia asteroides</i>		"	+++	+++	+++	+++	+++	+++	+++			
<i>Rhodococcus rhodochrous</i>		"	+++	+++	+++	+++	+++	+++	+++			
<i>Aerococcus viridans</i>	GT2116	Misc	+++	+++	+++	+++	+++	+++	+++			
<i>Fusobacterium necrophorum</i>	GT0238	Gm +	+++	+++	-	+++	+++	+++	+++			
<i>Fusobacterium prausnitzii</i>	ATCC27768	"	+++	-	+++	+++	+++	+++	+++			
<i>Gemella haemolysans</i>	GT2118	"	+++	+++	+++	+++	+++	+++	+++			
<i>Helicobacillus mobilis</i>		"	+++	++	+++	+++	+++	+++	+++			
<i>Phormidium ectocarpi</i>		Cyano	+	++	-	+++	+++	+++	+++			

TABLE 4 (cont'd): DOT BLOT HYBRIDIZATION OF 23S rRNA-TARGETED PROBES

Genus species	strain	div	PROBE HYBRIDIZATION				
			1730	1731	1658	1653	1654 1655 1651
Plectonema boryanum			++	++	-	+++	+++
Borrelia burgdorferi	Spiro		+++	++	++	-	+++
Borrelia turicatae	"		+++	+++	+++	-	+++
Leptospira interrogans-pomona	"		+++	+++	+++	-	+++
Leptospira biflexa (Patoc-Patoc)	"		+++	+++	+++	-	+++
Leptospira biflexa (CDC)	"		+++	+++	+++	-	+++
Spirochaeta aurantia	"		+++	+++	+++	-	+++
Bacteroides fragilis	25285 Bact		++	++	+	-	+++
Bacteroides fragilis	29771 "		++	++	+	-	+++
Bacteroides thetaiotaomicron	0572 "		+	+	+	-	+++
Bacteroides melaninogenicus	0011 "		-	-	+	++	+++
Bacteroides meningosepticum	0237 "		+	+	+	-	+++
Chlamydia psittaci	Chlam		-	-	-	-	-
Chlamydia trachomatis	LGV "		+++	+++	+	+++	+++
Chlorobium limicola	Misc.		+++	+++	+	++	+
Chloroflexus aurantiacus	Y400 "		+	+	+++	-	+
Deinococcus radiodurans	2608 "		+++	+++	+++	-	-
Planctomyces maris	2577 "		-	-	+++	-	-
Normal Stool RNA			+++	+++	+++	++	++
Mouse L-Cell			-	-	-	-	-
Wheat Germ			+	++	-	-	-
Normal Human Blood			-	-	-	-	-
Candida lusitanae			-	-	-	-	-
Candida parapsilosis			-	-	-	-	-
Candida tropicalis			-	-	-	-	-
Candida albicans			-	-	-	-	-
Candida albicans			-	-	-	-	-
Candida albicans			-	-	-	-	-

* Inclusivity and Exclusivity data was determined after overnight exposures.

** Each organism is represented by 100ng of CsTFA purified RNA. +++ = positive level of hybridization, + = barely detectable and - = zero, ND - not done.

TABLE 4 (cont'd): DOT BLOT HYBRIDIZATION OF 23S rRNA-TARGETED PROBES

Genus species	strain	div	PROBE HYBRIDIZATION									
			1512	1256	1398	1600	1601	1602	1596	1597		
Acinetobacter calcoaceticus	GT0002	Purple	++++	++++	+	++++	++++	++++	++++	++		
Aeromonas sobria	GT0007	gamma	++++	++++	++	++++	++++	++++	++++	++		
Alteromonas putrefaciens	GT1945	"	++++	++++	++	++++	++++	++++	++++	++		
Citrobacter amalonaticus	GT0690	"	++++	++++	+	++++	++++	++++	++++	++		
Citrobacter diversus	GT0030	"	++++	++++	+	++++	++++	++++	++++	++		
Citrobacter freundii	GT0687	"	++++	++++	+	++++	++++	++++	++++	++		
Edwardsiella tarda	GT0569	"	++++	++++	+	++++	++++	++++	++++	++		
Enterobacter agglomerans	GT0683	"	++++	++++	+	++++	++++	++++	++++	++		
Enterobacter cloacae	GT0686	"	++++	++++	+	++++	++++	++++	++++	++		
Enterobacter sakazakii	GT0062	"	++++	++++	++	++++	++++	++++	++++	++		
Escherichia coli	GT1665	"	++++	++++	+	++++	++++	++++	++++	++		
Escherichia coli	GT1592	"	++++	++++	+	++++	++++	++++	++++	++		
Escherichia coli	GT1659	"	++++	++++	+	++++	++++	++++	++++	++		
Haemophilus influenza	GT0244	"	++	+++	+	+++	+++	+++	+++	++		
Haemophilus influenza	ATCC33391	"	ND	ND	ND	ND	ND	ND	ND	++		
Haemophilus ducreyi	GT0243	"	++++	+++	++	++++	++++	++++	++++	++		
Hafnia alvei	GT0241	"	++++	+++	+	++++	++++	++++	++++	++		
Morganella morganii	GT0303	"	++++	+++	+	++++	++++	++++	++++	++		
Klebsiella pneumoniae	GT1500	"	++++	+++	+	++++	++++	++++	++++	++		
Proteus mirabilis	GT1496	"	++++	+++	+	++++	++++	++++	++++	++		
Providencia alcalifaciens	GT0371	"	++++	+++	+	++++	++++	++++	++++	++		
Pseudomonas aeruginosa	GT1909	"	++++	+++	++	++++	++++	++++	++++	++		
Salmonella arizona	GT0799	"	++++	+++	+	++++	++++	++++	++++	++		
Salmonella typhimurium	GT0389	"	++++	+++	+	++++	++++	++++	++++	++		
Serratia marcescens	GT0392	"	++++	+++	+	++++	++++	++++	++++	++		
Shigella flexneri	GT0798	"	++++	+++	++	++++	++++	++++	++++	++		
Vibrio parahaemolyticus	GT0568	"	++++	+++	++	++++	++++	++++	++++	++		
Xanthomonas maltophilia	GT0417	"	++++	+++	++	++++	++++	++++	++++	++		
Yersinia enterocolitica	GT0419	"	++++	+++	++	++++	++++	++++	++++	++		
Alcaligenes faecalis	GT0610	Purple	++++	+++	+	++++	++++	++++	++++	++		
Branhamella catarrhalis	GT0014	beta	++++	+++	+	++++	++++	++++	++++	++		
Chromobacterium violaceum	GT2022	"	++++	+++	+	++++	++++	++++	++++	++		
Kingella indologenes	0246	"	++++	+++	+	++++	++++	++++	++++	++		
Moraxella osloensis	GT0301	"	++++	+++	+	++++	++++	++++	++++	++		

TABLE 4 (cont'd): DOT BLOT HYBRIDIZATION OF 23S rRNA-TARGETED PROBES

Genus species	strain	div	PROBE HYBRIDIZATION									
			1512	1256	1398	1600	1601	1602	1596	1597		
<i>Morexella phenylpyruvica</i>	GT0302	"	+++	+++	+	+++	++	+++	+++	+++	+	+++
<i>Neisseria gonorrhoeae</i>	GT0315	"	+++	+++	+	+++	++	+++	+++	+++	+	+++
<i>Neisseria meningitidis</i>	GT0349	"	+++	+++	+	+++	+++	+++	+++	+++	+	+++
<i>Pseudomonas acidovorans</i>	GT0376	"	+++	+++	+	+++	+++	+++	+++	+++	+	+++
<i>Pseudomonas cepacia</i>	GT2015	"	+++	+++	+	+++	+++	+++	+++	+++	+	+++
<i>Rhodocyclus gelatinosa</i>	ATCC17013	"	+++	+++	+	+++	+++	+++	+++	+++	+	+++
<i>Vitreoscilla stercoraria</i>	GT0810	Purple	+++	+++	+	+++	+++	+++	+++	+++	+	+++
<i>Achromobacter xerosis</i>	ATCC33463	alpha	+++	+++	+	+++	+++	+++	+++	+++	+	+++
<i>Acidiphilium cryptum</i>	GT2021	"	+++	+++	+	+++	+++	+++	+++	+++	+	+++
<i>Agrobacterium tumefaciens</i>	ATCC23448	"	+++	+++	+	+++	+++	+++	+++	+++	+	+++
<i>Brucella abortus</i>	GT2025	"	+++	+++	+	+++	+++	+++	+++	+++	+	+++
<i>Flavobacterium capsulatum</i>	GT2023	"	+++	+++	+	+++	+++	+++	+++	+++	+	+++
<i>Mycoplana bullata</i>	GT2020	"	+++	+++	+	+++	+++	+++	+++	+++	+	+++
<i>Pseudomonas diminuta</i>	ATCC17023	"	+++	+++	+	+++	+++	+++	+++	+++	+	+++
<i>Rhodobacter sphaeroides</i>	ATCC25903	"	+++	+++	+	+++	+++	+++	+++	+++	+	+++
<i>Rhodospirillum rubrum</i>	ATCC25364	"	+++	+++	+	+++	+++	+++	+++	+++	+	+++
<i>Thiobacillus versutus</i>	ATCC 7757	delta	+++	+++	+	+++	+++	+++	+++	+++	+	+++
<i>Desulfovibrio desulfuricans</i>	GT2095	Gm -	+++	+++	+	+++	+++	+++	+++	+++	+	+++
<i>Cardiobacterium hominis</i>	GT2172	Gm -	+++	+++	+	+++	+++	+++	+++	+++	+	+++
<i>Francisella tularensis</i>	GT0022	Campy	+++	+++	+	+++	+++	+++	+++	+++	+	+++
<i>Campylobacter jejuni</i>	GT0026	"	+++	+++	+	+++	+++	+++	+++	+++	+	+++
<i>Campylobacter pylori</i>	GT0027	"	+++	+++	+	+++	+++	+++	+++	+++	+	+++
<i>Campylobacter sputorum</i>	GT0803	lowG+C	+++	+++	+	+++	+++	+++	+++	+++	+	+++
<i>Bacillus brevis</i>	GT0804	Gm +	+++	+++	+	+++	+++	+++	+++	+++	+	+++
<i>Bacillus subtilis</i>	ATCC25537	"	+++	+++	+	+++	+++	+++	+++	+++	+	+++
<i>Clostridium clostridioforme</i>	GT0567	"	+++	+++	+	+++	+++	+++	+++	+++	+	+++
<i>Clostridium sordellii</i>	ATCC 3587	"	+++	+++	+	+++	+++	+++	+++	+++	+	+++
<i>Clostridium sporogenes</i>	ATCC19401	"	+++	+++	+	+++	+++	+++	+++	+++	+	+++
<i>Clostridium histolyticum</i>	ATCC13124	"	+++	+++	+	+++	+++	+++	+++	+++	+	+++
<i>Clostridium perfringens</i>	ATCC2582	"	+++	+++	+	+++	+++	+++	+++	+++	+	+++
<i>Clostridium ramosum</i>	ATCC33403	"	+++	+++	+	+++	+++	+++	+++	+++	+	+++
<i>Kurthia zopfii</i>	GT0256	"	+++	+++	+	+++	+++	+++	+++	+++	+	+++
<i>Lactobacillus acidophilus</i>	GT0258	"	+++	+++	+	+++	+++	+++	+++	+++	+	+++
<i>Lactobacillus plantarum</i>		"	+++	+++	+	+++	+++	+++	+++	+++	+	+++

TABLE 4 (cont'd): DOT BLOT HYBRIDIZATION OF 23S rRNA-TARGETED PROBES

Genus species	strain	div	PROBE HYBRIDIZATION											
			1512	1256	1398	1600	1601	1602	1596	1597				
<i>Listeria monocytogenes</i>	IG3299	"	+++	+++	+	+++	+++	+++	++	+++				
<i>Mycoplasma pneumoniae</i>	ATCC1531	"	+++	-	+	++	-	-	-	-				
<i>Mycoplasma putrefaciens</i>	ATCC15718	"	+++	-	-	++	++	+++	+++	+++				
<i>Peptostreptococcus productus</i>	ATCC27340	"	+++	+++	-	++	++	+++	+++	+++				
<i>Planococcus citreus</i>	ATCC14404	"	+++	+++	+	+++	+++	+++	+++	+++				
<i>Staphylococcus aureus</i>	GT0399	"	+++	+++	+	+++	+++	+++	+++	+++				
<i>Staphylococcus aureus</i>	GT1711	"	+++	+++	+	+++	+++	+++	+++	+++				
<i>Staphylococcus epidermidis</i>	GT0401	"	+++	+++	+	+++	+++	+++	+++	+++				
<i>Staphylococcus haemolyticus</i>	ATCC29970	"	+++	+++	+	+++	+++	+++	+++	+++				
<i>Streptococcus agalactiae</i>	GT0405	"	+++	+++	+	+++	+++	+++	+++	+++				
<i>Streptococcus bovis</i>	GT0658	"	+++	+++	+	+++	+++	+++	+++	+++				
<i>Streptococcus faecalis</i>	GT0406	"	+++	+++	+	+++	+++	+++	+++	+++				
<i>Streptococcus morbillorum</i>	GT2194	"	+++	+++	+	+++	+++	+++	+++	+++				
<i>Streptococcus mutans</i>	GT0412	"	+++	+++	+	+++	+++	+++	+++	+++				
<i>Streptococcus pneumoniae</i>	GT0408	"	+++	+++	+	+++	+++	+++	+++	+++				
<i>Streptococcus salivarius</i>	GT0410	"	+++	+++	+	+++	+++	+++	+++	+++				
<i>Streptococcus sanguis</i>	GT0411	"	ND	ND	+	+++	+++	+++	+++	+++				
<i>Bifidobacterium dentium</i>	GT0012	hiG+C	-	+++	+	+++	+++	+++	+++	+++				
<i>Corynebacterium genitalium</i>	GT0045	Gm +	++	+++	+	+++	+++	+++	+++	+++				
<i>Corynebacterium glutamicum</i>	GT2120	"	++	+++	+	+++	+++	+++	+++	+++				
<i>Corynebacterium pseudodiphtheriticum</i>	GT2119	"	+++	+++	+	+++	+++	+++	+++	+++				
<i>Corynebacterium pseudotuberculosis</i>	GT2122	"	+++	+++	+	+++	+++	+++	+++	+++				
<i>Corynebacterium pyogenes</i>	GT2121	"	++	+++	+	+++	+++	+++	+++	+++				
<i>Corynebacterium xerosis</i>	GT0046	"	+	++	+	+++	+++	+++	+++	+++				
<i>Corynebacterium xerosis</i>	BCG	"	+	++	+	+++	+++	+++	+++	+++				
<i>Mycobacterium bovis</i>	GT2191	"	+++	+++	+	+++	+++	+++	+++	+++				
<i>Mycobacterium kansasii</i>	GT2191	"	+++	+++	+	+++	+++	+++	+++	+++				
<i>Nocardia asteroides</i>	GT2116	Misc	++	+++	+	+++	+++	+++	+++	+++				
<i>Rhodococcus rhodochrous</i>	GT0238	Gm +	++	+++	+	+++	+++	+++	+++	+++				
<i>Aerococcus viridans</i>	ATCC27758	"	+++	+++	+	+++	+++	+++	+++	+++				
<i>Fusobacterium necrophorum</i>	GT2118	"	+++	+++	+	+++	+++	+++	+++	+++				
<i>Fusobacterium prausnitzii</i>	GT2118	"	+++	+++	+	+++	+++	+++	+++	+++				
<i>Gemella haemolysans</i>		"	+++	+++	+	+++	+++	+++	+++	+++				
<i>Hellobacillus mobilis</i>		"	+++	+++	+	+++	+++	+++	+++	+++				
<i>Phormidium ectocarpi</i>		Cyano	+++	+++	+	+++	+++	+++	+++	+++				

TABLE 4 (cont'd): DOT BLOT HYBRIDIZATION OF 23S rRNA-TARGETED PROBES

Genus species	strain	div	PROBE HYBRIDIZATION									
			1512	1256	1398	1600	1601	1602	1596	1597		
Plectonema boryanum			+++	+++	+	+	+++	+++	++	++		
Borrelia burgdorferi	Spiro		+++	+++	+	+++	+++	+++	+++	++		
Borrelia turicatae	"		+++	+++	+	+++	+++	+++	+++	++		
Leptospira interrogans-pomona	"		+++	+++	+	+++	+++	+++	+++	++		
Leptospira biflexa (Patoc-Patoc)	"		+++	+++	+	+++	+++	+++	+++	++		
Leptospira biflexa (CDC)	"		+++	+++	+	+++	+++	+++	+++	++		
Spirochaeta aurantia			+++	+++	+	+++	+++	+++	+++	++		
Bacteroides fragilis	25285 Bact		+++	+++	+	+++	+++	+++	+++	++		
Bacteroides fragilis	29771 "		+++	+++	+	+++	+++	+++	+++	++		
Bacteroides thetaiotaomicron	0572 "		+++	+++	+	+	+++	+++	+++	++		
Bacteroides melaninogenicus	0011 "		+++	+++	+	+	+++	+++	+++	++		
Bacteroides meningosepticum	0237 "		+++	+++	+	+++	+++	+++	+++	++		
Flavobacterium meningosepticum			+	+++	-	+++	+++	+++	+++	++		
Chlamydia psittaci	LGV		+	+++	-	+++	+++	+++	+++	++		
Chlamydia trachomatis	Misc.		+++	+++	+	+++	+++	+++	+++	++		
Chlorobium limicola			+	+++	+	+	+++	+++	+++	++		
Chloroflexus aurantiacus	Y400		+++	+++	+	+++	+++	+++	+++	++		
Deinococcus radiodurans	2608 "		+++	+++	+	+++	+++	+++	+++	++		
Planctomyces maris	2577		+++	+++	+	+++	+++	+++	+++	++		
Normal Stool RNA			+++	+++	-	+++	+++	+++	+++	++		
Mouse L-Cell			-	-	-	-	-	-	-	-		
Wheat Germ			-	-	-	-	-	-	-	-		
Normal Human Blood			+	-	-	-	-	-	-	-		
Candida lusitanae	403-87		+	-	-	-	-	-	-	-		
Candida parapsilosis	882-88		+	-	-	-	-	-	-	-		
Candida tropicalis	224-87		+	-	-	-	-	-	-	-		
Candida albicans	1008-88		+	-	-	-	-	-	-	-		
Candida albicans	223-87		+	-	-	-	-	-	-	-		
Candida albicans	819-88		+	-	-	-	-	-	-	-		

* Inclusivity and Exclusivity data was determined after overnight exposures.

** Each organism is represented by 100ng of CsTFA purified RNA. +++ = positive level of hybridization, + = barely detectable and - = zero, ND = not done.

TABLE 5: DOT BLOT HYBRIDIZATION - GRAM POSITIVE & GRAM NEGATIVE PROBES

Genus species	strain	div	PROBE HYBRIDIZATION						
			16S RNA-TARGET		23S RNA-TARGET				
			1744	1745	1746	1657	1656	1598	1599 1595
Acinetobacter calcoaceticus	GT0002	Purple	-	-	-	+++	-	++	+++
Aeromonas sobria	GT0007	gamma	-	-	-	+++	-	++	+++
Alteromonas putrefaciens	GT1945	"	-	-	-	+++	-	+++	+++
Citrobacter amalonaticus	GT0690	"	-	-	-	+++	-	+++	+++
Citrobacter diversus	GT0030	"	-	-	-	+++	-	+++	+++
Citrobacter freundii	GT0687	"	-	-	-	+++	-	+++	+++
Edwardsiella tarda	GT0569	"	-	-	-	+++	-	+++	+++
Enterobacter agglomerans	GT0683	"	-	-	-	+++	-	++	+++
Enterobacter cloacae	GT0686	"	+	-	-	+++	-	++	+++
Enterobacter sakazakii	GT0062	"	-	-	-	+++	-	+++	+++
Escherichia coli	GT1665	"	-	-	-	+++	-	+++	+++
Escherichia coli	GT1592	"	-	-	-	+++	-	+++	+++
Echerichia coli	GT1659	"	-	-	-	+++	+	++	+++
Haemophilus influenza	ATCC33391	"	-	-	-	+++	+++	+++	+++
Haemophilus ducreyi	GT0243	"	-	-	-	+++	-	+++	+++
Haemophilus alvei	GT0241	"	-	-	-	+++	-	+++	+++
Morganella morganii	GT0303	"	-	-	-	+++	-	+++	+++
Klebsiella pneumoniae	GT1500	"	-	-	-	+++	-	+++	+++
Proteus mirabilis	GT1496	"	-	-	-	+++	-	+++	+++
Providencia alcalifaciens	GT0371	"	-	-	-	+++	-	++	+++
Pseudomonas aeruginosa	GT1909	"	-	-	-	+++	-	++	+++
Salmonella arizona	GT0799	"	-	-	-	+++	-	++	+++
Salmonella typhimurium	GT0389	"	-	-	-	+++	-	+++	+++
Serratia marcescens	GT0392	"	-	-	-	+++	-	+++	+++
Shigella flexneri	GT0798	"	-	-	-	+++	-	+++	+++
Vibrio parahaemolyticus	GT0568	"	-	-	-	+++	-	+++	+++
Xanthomonas maltophilia	GT0417	"	-	-	-	+++	+	+++	+++
Yersinia enterocolitica	GT0419	"	-	-	-	+++	-	+++	+++
Alcaligenes faecalis	GT0610	Purple	-	-	-	+++	+++	+++	+++
Branhamella catarrhalis	GT0014	beta	-	-	-	+++	-	+++	+++
Chromobacterium violaceum	GT2022	"	-	-	-	+++	++	+++	+++
Kingella indologenes	0246	"	-	-	-	+++	-	+++	+++
Moraxella osloensis	GT0301	"	-	-	-	+++	-	+++	+++
Morexella phenylpyruvica	GT0302	"	-	-	-	+++	-	+++	+++

TABLE 5: (cont'd) DOT BLOT HYBRIDIZATION - GRAM POSITIVE & GRAM NEGATIVE PROBES

Genus species	strain	div	PROBE HYBRIDIZATION							
			16S RNA-TARGET				23S RNA-TARGET			
			1744	1745	1746	1657	1656	1598	1599	1595
Borrelia burgdorferi		Spiro	-	+	+++	-	++++	++	-	+
Borrelia turicatae		"	-	++	++++	-	++++	++++	-	+
Leptospira interrogans-pomona		"	-	+	-	-	++++	++++	++	+
Leptospira biflexa (Patoc-Patoc)		"	-	-	-	-	++++	++++	++	+
Leptospira biflexa (CDC)		"	-	-	-	-	++++	++++	++	+
Spirochaeta aurantia		"	++++	-	-	-	++++	++	++++	+
Bacteroides fragilis	25285 Bact	"	++++	-	-	-	-	-	++++	+
Bacteroides fragilis	29771 "	"	++++	-	-	-	-	-	++++	+
Bacteroides thetaiotaomicron	0572 "	"	+	-	-	-	-	-	++++	+
Bacteroides melaninogenicus	0011 "	"	++++	-	-	-	-	-	++++	+
Flavobacterium meningosepticum	0237 "	"	++++	-	-	-	+	-	++++	-
Chlamydia psittaci		Chlamy	-	-	-	-	-	+	++	++
Chlamydia trachomatis	LGV "	"	-	-	-	-	-	++	-	++
Chlorobium limicola		Misc	++++	-	-	-	+++	-	-	++
Chloroflexus aurantiacus	Y400 "	"	++++	-	-	-	+++	-	++++	+
Deinococcus radiodurans	2608 "	"	+++	-	-	-	+++	-	++++	+
Planctomyces maris	2577 "	"	-	-	-	-	-	-	-	-
Normal Stool RNA			+++	+++	++	++	+++	+	+++	++
Mouse L-Cell			-	-	-	-	-	-	-	-
Wheat Germ			-	-	-	-	-	-	-	-
Normal Human Blood			-	-	-	-	-	-	-	-
Candida lusitanae	403-87		-	-	-	-	-	-	-	-
Candida parapsilosis	882-88		-	-	-	-	-	-	-	-
Candida tropicalis	224-87		-	-	-	-	-	-	-	-
Candida albicans	1008-88		-	-	-	-	-	-	-	-
Candida albicans	223-87		-	-	-	-	-	-	-	-
Candida albicans	819-88		-	-	-	-	-	-	-	-

* Inclusivity and Exclusivity data was determined after overnight exposures.

** Each organism is represented by 100ng of CstFA purified RNA.

*** +++ = positive level of hybridization, + = barely detectable and - = zero. Probe 1746 was hybridized at 37 C and washes at 50 C.

TABLE 5: (cont'd) DOT BLOT HYBRIDIZATION - GRAM POSITIVE & GRAM NEGATIVE PROBES

Genus species	strain	div	PROBE HYBRIDIZATION						
			16S RNA-TARGET		23S RNA-TARGET				
			1744	1745	1746	1657	1656	1598	1599 1595
Mycoplasma pneumoniae	ATCC15531	"	+++	+++	+++	-	-	+++	+++
Mycoplasma putrefaciens	ATCC15718	"	+++	+++	+++	-	-	+++	+++
Peptostreptococcus productus	ATCC27340	"	+++	+++	+++	-	+	+++	+++
Planococcus citreus	ATCC14404	"	+++	+++	+++	-	+	+++	+++
Staphylococcus aureus	GT0399	"	+++	+++	+++	-	+	+++	+++
Staphylococcus aureus	GT1711	"	+++	+++	+++	-	+	+++	+++
Staphylococcus epidermidis	GT0401	"	+++	+++	+++	-	+	+++	+++
Staphylococcus haemolyticus	ATCC29970	"	+++	+++	+++	-	+	+++	+++
Streptococcus agalactiae	GT0405	"	+++	+++	+++	-	+	+++	+++
Streptococcus bovis	GT0668	"	+++	+++	+++	-	+	+++	+++
Streptococcus faecalis	GT0406	"	+++	+++	+++	-	+	+++	+++
Streptococcus morbillorum	GT2194	"	+++	+++	+++	-	+	+++	+++
Streptococcus mutans	GT0412	"	+++	+++	+++	-	+	+++	+++
Streptococcus pneumoniae	GT0408	"	+++	+++	+++	-	+	+++	+++
Streptococcus salivarius	GT0410	"	+++	+++	+++	-	+	+++	+++
Streptococcus sanguis	GT0411	"	+++	+++	+++	-	+	+++	+++
Streptococcus dentium	GT0012	h1G+C	+++	+++	+++	-	+	+++	+++
Bifidobacterium dentium	GT0045	Gm +	+++	+++	+++	-	+	+++	+++
Corynebacterium genitalium	GT2120	"	+++	+++	+++	-	+	+++	+++
Corynebacterium glutamicum	GT2119	"	+++	+++	+++	-	+	+++	+++
Corynebacterium pseudodiphtheriticum	GT2122	"	+++	+++	+++	-	+	+++	+++
Corynebacterium pseudotuberculosis	GT2121	"	+++	+++	+++	-	+	+++	+++
Corynebacterium pyogenes	GT0046	"	+++	+++	+++	-	+	+++	+++
Corynebacterium xerosis	BCG	"	+++	+++	+++	-	+	+++	+++
Corynebacterium bovis	GT2191	"	+++	+++	+++	-	+	+++	+++
Mycobacterium kansasii	GT2116	Misc	+++	+++	+++	-	+	+++	+++
Nocardia asteroides	GT0238	Gm +	+++	+++	+++	-	+	+++	+++
Rhodococcus rhodochrous	ATCC27768	"	+++	+++	+++	-	+	+++	+++
Aerococcus viridans	GT2118	"	+++	+++	+++	-	+	+++	+++
Fusobacterium necrophorum	ATCC27768	"	+++	+++	+++	-	+	+++	+++
Fusobacterium prausnitzii	GT2118	"	+++	+++	+++	-	+	+++	+++
Gemella haemolysans	"	"	+++	+++	+++	-	+	+++	+++
Helicobacter mobilis	"	"	+++	+++	+++	-	+	+++	+++
Phormidium ectocarpi	"	"	+++	+++	+++	-	+	+++	+++
Plectonema boryanum	"	Cyano	+++	+++	+++	-	+	+++	+++

TABLE 5: (cont'd) DOT BLOT HYBRIDIZATION - GRAM POSITIVE & GRAM NEGATIVE PROBES

Genus species	strain	div	PROBE HYBRIDIZATION							
			16S RNA-TARGET				23S RNA-TARGET			
			1744	1745	1746	1657	1656	1598	1599	1595
<i>Neisseria gonorrhoeae</i>	GT0315	"	+	+	+	+	+	+	+	+
<i>Neisseria meningitidis</i>	GT0349	"	+	+	+	+	+	+	+	+
<i>Pseudomonas acidovorans</i>	GT0376	"	+	+	+	+	+	+	+	+
<i>Pseudomonas cepacia</i>	GT2015	"	+	+	+	+	+	+	+	+
<i>Rhodocyclus gelatinosa</i>	ATCC17013	"	+	+	+	+	+	+	+	+
<i>Vitreoscilla stercoraria</i>	GT0810	Purple	+	+	+	+	+	+	+	+
<i>Achromobacter xerosis</i>	ATCC33463	alpha	+	+	+	+	+	+	+	+
<i>Acidiphilium cryptum</i>	GT2021	"	+	+	+	+	+	+	+	+
<i>Agrobacterium tumefaciens</i>	ATCC23448	"	+	+	+	+	+	+	+	+
<i>Brucella abortus</i>	GT2025	"	+	+	+	+	+	+	+	+
<i>Flavobacterium capsulatum</i>	GT2023	"	+	+	+	+	+	+	+	+
<i>Mycoplasma bullata</i>	GT2020	"	+	+	+	+	+	+	+	+
<i>Pseudomonas diminuta</i>	ATCC17023	"	+	+	+	+	+	+	+	+
<i>Rhodobacter sphaeroides</i>	ATCC25903	"	+	+	+	+	+	+	+	+
<i>Rhodospirillum rubrum</i>	ATCC25364	"	+	+	+	+	+	+	+	+
<i>Thiobacillus versutus</i>	ATCC 7757	delta	+	+	+	+	+	+	+	+
<i>Desulfovibrio desulfuricans</i>	GT2095	Misc	+	+	+	+	+	+	+	+
<i>Cardiobacterium hominis</i>	GT2172	Gm -	+	+	+	+	+	+	+	+
<i>Francisella tularensis</i>	GT0022	Campy	+	+	+	+	+	+	+	+
<i>Campylobacter jejuni</i>	GT0026	"	+	+	+	+	+	+	+	+
<i>Campylobacter pylori</i>	GT0027	"	+	+	+	+	+	+	+	+
<i>Campylobacter sputorum</i>	GT0803	lowG+C	+	+	+	+	+	+	+	+
<i>Bacillus brevis</i>	GT0804	Gm +	+	+	+	+	+	+	+	+
<i>Bacillus subtilis</i>	ATCC25537	"	+	+	+	+	+	+	+	+
<i>Clostridium clostridioforme</i>	GT0567	"	+	+	+	+	+	+	+	+
<i>Clostridium sordellii</i>	ATCC 3587	"	+	+	+	+	+	+	+	+
<i>Clostridium sporogenes</i>	ATCC19401	"	+	+	+	+	+	+	+	+
<i>Clostridium histolyticum</i>	ATCC13124	"	+	+	+	+	+	+	+	+
<i>Clostridium perfringens</i>	ATCC25582	"	+	+	+	+	+	+	+	+
<i>Clostridium ramosum</i>	ATCC33403	"	+	+	+	+	+	+	+	+
<i>Kurthia zopfii</i>	GT0256	"	+	+	+	+	+	+	+	+
<i>Lactobacillus acidophilus</i>	GT0258	"	+	+	+	+	+	+	+	+
<i>Lactobacillus plantarum</i>	IC3299	"	+	+	+	+	+	+	+	+
<i>Listeria monocytogenes</i>			+	+	+	+	+	+	+	+

Claims

1. A nucleic acid sequence capable of hybridising to rRNA or rDNA of eubacteria, but not to rRNA or rDNA of Mouse L cells, wheat germ, human blood or *Candida albicans*, the sequence being complementary to, or homologous to, at least 90% of a sequence comprising any ten consecutive nucleotides within a probe which is:

- Probe 1661: 5'-TATTACCGCGGCTGCTGGCACGGAGTTAGCCG-3' ;
 Probe 1739: 5'-GCGTGGACTACCGGGGTATCTAATCCTGTTTGCTCCCCACGCTTTCG-3' ;
 Probe 1746: 5'-TCATAAGGGGCATGATGATTTGACGTCAT-3' ;
 Probe 1743: 5'-GTACAAGGCCCGGGAACGTATTCACCG-3' ;
 5 Probe 1639: 5'-ACGGTTACCTTGTTACGACTT-3' ;
 Probe 1640: 5'-ACGGCTACCTTGTTACGACTT-3' ;
 Probe 1641: 5'-ACGGATACCTTGTTACGACTT-3' ;
 Probe 1656: 5'-CTACCTGTGTCGGTTTGCGGTACGGGC-3' ;
 Probe 1657: 5'-GGTATTCTCTACCTGACCACCTGTGTCGGTTTGGGGTACG-3' ;
 10 Probe 1653: 5'-CCTTCTCCCGAAGTTACGGGGGCATTTTGCCTAGTTCCTT-3' ;
 Probe 1654: 5'-CCTTCTCCCGAAGTTACGGGGTCATTTTGCCGAGTTCCTT-3' ;
 Probe 1655: 5'-CCTTCTCCCGAAGTTACGGCACCATTGTCGGAGTTCCTT-3' ;
 Probe 1651: 5'-CTCCTCTTAACCTTCCAGCACCGGGCAGGC-3' ;
 Probe 1652: 5'-TTCGATCAGGGGCTTCGCTTGCGCTGACCCCATCAATTAA-3' ;
 15 Probe 1595: 5'-CGATATGAACCTTTGGGCGGTATCAGCCTGTTATCCCCGG-3' ;
 Probe 1600: 5'-CAGCCCCAGGATGAGATGAGCCGACATCGAGGTGCCAAAC-3' ;
 Probe 1601: 5'-CAGCCCCAGGATGTGATGAGCCGACATCGAGGTGCCAAAC-3' ;
 Probe 1602: 5'-CAGCCCCAGGATGCGATGAGCCGACATCGAGGTGCCAAAC-3' ;
 Probe 1598: 5'-CGTACCGCTTTAAATGGCGAACAGCCATACCCTTGGGACC-3' ;
 20 Probe 1599 : 5'-CGTGCCGCTTTAATGGGCGAACAGCCCAACCCTTGGGACC-3' ;
 Probe 1596: 5'-GATAGGGACCGAACTGTCTCACGACGTTTGAACCCAGCT-3' ; or
 Probe 1597: 5'-GATAGGGACCGAACTGTCTCACGACGTTCTGAACCCAGCT-3' .
2. A nucleic acid sequence as claimed in claim 1 which is probe 1661, probe 1739, probe 1746, probe
 25 1743, probe 1639, probe 1640, probe 1641, probe 1656, probe 1657, probe 1653, probe 1654, probe
 1655, probe 1651, probe 1652, probe 1595, probe 1600, probe 1601, probe 1602, probe 1598, probe
 1599, probe 1596 or probe 1597, or a complementary sequence thereto.
3. A set of probes comprising at least two nucleic acid sequences at least one of which is 1661, 1739,
 30 1746, 1743, 1639, 1640, 1641, 1656, 1657, 1653, 1654, 1655, 1651, 1652, 1595, 1600, 1601, 1602,
 1598, 1599, 1596 or 1597, or a complementary sequence thereto.
4. A method of detecting the presence of eubacteria in a sample, the method comprising:
 35 (a) contacting a sample with at least one nucleic acid sequence which is probe 1661, 1739, 1746,
 1743, 1639, 1640, 1641, 1656, 1657, 1653, 1654, 1655, 1651, 1652, 1595, 1600, 1601, 1602, 1598,
 1599, 1596 or 1597, or a complementary sequence thereto, under conditions that allow the
 sequence to hybridise to rRNA or rDNA of the eubacteria, if present in the sample, but not to form
 hybrid nucleic acid complexes to rRNA or rDNA of non-eubacteria; and
 40 (b) detecting the hybrid nucleic acid complexes as an indication of the presence of eubacteria in the
 sample.
5. A method as claimed in claim 4 wherein the eubacteria is gram-positive or gram-negative and the
 nucleic acid sequence is the probe 1599, 1656 or 1746.

45 Patentansprüche

1. Nucleinsäuresequenz, die mit rRNA oder rDNA von Eubakterien hybridisieren kann, aber nicht mit
 rRNA oder rDNA von Maus-L-Zellen, Weizenkeimen, menschlichem Blut oder Candida albicans, wobei
 50 die Sequenz komplementär zu oder homolog mit mindestens 90 % einer Sequenz ist, die beliebige
 zehn aufeinanderfolgende Nucleotide innerhalb einer Sonde umfaßt, nämlich:
 Sonde 1661: 5'-TATTACCGCGGCTGCTGGCACGGAGTTAGCCG-3' ;
 Sonde 1739: 5'-GCGTGGACTACCGGGGTATCTAATCCTGTTTGCTCCCCACGCTTTCG-3' ;
 Sonde 1746: 5'-TCATAAGGGGCATGATGATTTGACGTCAT-3' ;
 Sonde 1743: 5'-GTACAAGGCCCGGGAACGTATTCACCG-3' ;
 55 Sonde 1639: 5'-ACGGTTACCTTGTTACGACTT-3' ;
 Sonde 1640: 5'-ACGGCTACCTTGTTACGACTT-3' ;
 Sonde 1641: 5'-ACGGATACCTTGTTACGACTT-3' ;
 Sonde 1656: 5'-CTACCTGTGTCGGTTTGCGGTACGGGC-3' ;

Sonde 1657: 5'-GGTATTCTCTACCTGACCACCTGTGTGCGGTTTGGGGTACG-3';
 Sonde 1653: 5'-CCTTCTCCCGAAGTTACGGGGGCATTTTGCCTAGTTCCTT-3';
 Sonde 1654: 5'-CCTTCTCCCGAAGTTACGGGGGCATTTTGCCGAGTTCCTT-3';
 Sonde 1655: 5'-CCTTCTCCCGAAGTTACGGCACCATTGCGGAGTTCCTT-3';
 5 Sonde 1651: 5'-CTCCTCTTAACCTTCCAGCACCGGGCAGGC-3';
 Sonde 1652: 5'-TTCGATCAGGGGCTTCGCTTGCCTGACCCCATCAATTAA-3';
 Sonde 1595: 5'-CGATATGAACTCTGGGCGGTATCAGCCTGTTATCCCCGG-3';
 Sonde 1600: 5'-CAGCCCCAGGATGAGATGAGCCGACATCGAGGTGCCAAAC-3';
 Sonde 1601: 5'-CAGCCCCAGGATGTGATGAGCCGACATCGAGGTGCCAAAC-3';
 10 Sonde 1602: 5'-CAGCCCCAGGATGCGATGAGCCGACATCGAGGTGCCAAAC-3';
 Sonde 1598: 5'-CGTACCGCTTTAAATGGCGAACAGCCATACCCTTGGGACC-3';
 Sonde 1599: 5'-CGTACCGCTTTAAATGGCGAACAGCCCAACCCTTGGGACC-3';
 Sonde 1596: 5'-GATAGGGACCGAACTGTCTCAGACCTTTGAACCCAGCT-3'; oder
 Sonde 1597: 5'-GATAGGGACCGAACTGTCTCAGACGTTCTGAACCCAGCT-3'.

- 15 2. Nucleinsäuresequenz nach Anspruch 1, nämlich Sonde 1661, Sonde 1739, Sonde 1746, Sonde 1743, Sonde 1639, Sonde 1640, Sonde 1641, Sonde 1656, Sonde 1657, Sonde 1653, Sonde 1654, Sonde 1655, Sonde 1651, Sonde 1652, Sonde 1595, Sonde 1600, Sonde 1601, Sonde 1602, Sonde 1598, Sonde 1599, Sonde 1596 oder Sonde 1597 oder eine komplementäre Sequenz dazu.
- 20 3. Sondensatz, umfassend mindestens zwei Nucleinsäuresequenzen, wobei mindestens eine dieser 1661, 1739, 1746, 1743, 1639, 1640, 1641, 1656, 1657, 1653, 1654, 1655, 1651, 1652, 1595, 1600, 1601, 1602, 1598, 1599, 1596 oder 1597 oder eine komplementäre Sequenz dazu ist.
- 25 4. Verfahren zum Nachweis der Gegenwart von Eubakterien in einer Probe, umfassend:
 (a) Inkontaktbringen einer Probe mit mindestens einer Nucleinsäuresequenz, die Sonde 1661, 1739, 1746, 1743, 1639, 1640, 1641, 1656, 1657, 1653, 1654, 1655, 1651, 1652, 1595, 1600, 1601, 1602, 1598, 1599, 1596 oder 1597 oder eine komplementäre Sequenz dazu ist, unter Bedingungen, bei
 30 denen die Sequenz mit rRNA oder rDNA von Eubakterien, falls sie in der Probe vorhanden sind, hybridisieren, aber keine Hybridnucleinsäurekomplexe mit rRNA oder rDNA von Nicht-Eubakterien bilden kann; und
 (b) Nachweis der Hybridnucleinsäurekomplexe als ein Anzeichen für die Gegenwart von Eubakterien in der Probe.
- 35 5. Verfahren nach Anspruch 4, wobei die Eubakterien gram-positive oder gram-negative Bakterien sind und die Nucleinsäuresequenz die Sonde 1599, 1656 oder 1746 ist.

Revendications

- 40 1. Séquence d'acide nucléique capable de s'hybrider à l'ARNr ou l'ADNr d'eubactéries, mais non à l'ARNr ou l'ADNr de cellules L de souris, de germe de blé, de sang humain ou de Candida albicans, la séquence étant complémentaire de, ou homologue à, au moins 90% d'une séquence comprenant n'importe quelle dizaine de nucléotides consécutifs dans une sonde qui est :
- 45 Sonde 1661: 5'-TATTACCGCGGCTGCTGGCACGGAGTTAGCCG-3';
 Sonde 1739: 5'-GCGTGGACTACCGGGGTATCTAATCCTGTTTGCTCCCCACGCTTTCG-3';
 Sonde 1746: 5'-TCATAAGGGGCATGATGATTTGACGTCAT-3';
 Sonde 1743: 5'-GTACAAGGCCCGGGAACGTATTCACCG-3';
 Sonde 1639: 5'-ACGGTTACCTTGTTACGACTT-3';
 Sonde 1640: 5'-ACGGCTACCTTGTTACGACTT-3';
 50 Sonde 1641: 5'-ACGGATACCTTGTTACGACTT-3';
 Sonde 1656: 5'-CTACCTGTGTGCGGTTTGCGGTAACGGGC-3';
 Sonde 1657: 5'-GGTATTCTCTACCTGACCACCTGTGTGCGGTTTGGGGTACG-3';
 Sonde 1653: 5'-CCTTCTCCCGAAGTTACGGGGGCATTTTGCCTAGTTCCTT-3';
 Sonde 1654: 5'-CCTTCTCCCGAAGTTACGGGGGCATTTTGCCGAGTTCCTT-3';
 55 Sonde 1655: 5'-CCTTCTCCCGAAGTTACGGCACCATTGCGGAGTTCCTT-3';
 Sonde 1651: 5'-CTCCTCTTAACCTTCCAGCACCGGGCAGGC-3';
 Sonde 1652: 5'-TTCGATCAGGGGCTTCGCTTGCCTGACCCCATCAATTAA-3';
 Sonde 1595: 5'-CGATATGAACTCTGGGCGGTATCAGCCTGTTATCCCCGG-3';

Sonde 1600: 5'-CAGCCCCAGGATGAGATGAGCCGACATCGAGGTGCCAAAC-3' ;
 Sonde 1601: 5'-CAGCCCCAGGATGTGATGAGCCGACATCGAGGTGCCAAAC-3' ;
 Sonde 1602: 5'-CAGCCCCAGGATGCGATGAGCCGACATCGAGGTGCCAAAC-3' ;
 Sonde 1598: 5'-CGTACCGCTTTAAATGGCGAACAGCCATACCCTTGGGACC-3' ;
 Sonde 1599: 5'-CGTGCCGCTTTAATGGGCGAACAGCCCAACCCTTGGGACC-3' ;
 Sonde 1596: 5'-GATAGGGACCGAACTGTCTCAGGACGTTTTGAACCCAGCT-3' ; ou
 Sonde 1597: 5'-GATAGGGACCGAACTGTCTCAGGTTCTGAACCCAGCT-3' .

2. Séquence d'acide nucléique selon la revendication 1 qui est la sonde 1661, la sonde 1739, la sonde 1746, la sonde 1743, la sonde 1639, la sonde 1640, la sonde 1641, la sonde 1656, la sonde 1657, la sonde 1653, la sonde 1654, la sonde 1655, la sonde 1651, la sonde 1652, la sonde 1595, la sonde 1600, la sonde 1601, la sonde 1602, la sonde 1598, la sonde 1599, la sonde 1596 ou la sonde 1597, ou une séquence complémentaire de l'une de ces sondes.
3. Jeu de sondes comprenant au moins deux séquences d'acide nucléique dont l'une au moins est 1661, 1739, 1746, 1743, 1639, 1640, 1641, 1656, 1657, 1653, 1654, 1655, 1651, 1652, 1595, 1600, 1601, 1602, 1598, 1599, 1596 ou 1597, ou une séquence complémentaire de celle-ci.
4. Procédé de détection de la présence d'eubactéries dans un échantillon, qui comprend :
 - (a) la mise en contact d'un échantillon avec au moins une séquence d'acide nucléique qui est la sonde 1661, 1739, 1746, 1743, 1639, 1640, 1641, 1656, 1657, 1653, 1654, 1655, 1651, 1652, 1595, 1600, 1601, 1602, 1598, 1599, 1596 ou 1597, ou une séquence complémentaire de celle-ci, dans des conditions qui permettent à la séquence de s'hybrider à l'ARNr ou l'ADNr des eubactéries, si elles sont présentes dans l'échantillon mais qui ne permettent pas de former des complexes d'acides nucléiques hybrides avec l'ARNr ou l'ADNr de non-eubactéries; et
 - (b) la détection des complexes d'acides nucléiques hybrides en tant qu'indication de la présence d'eubactéries dans l'échantillon.
5. Procédé selon la revendication 4, dans lequel les eubactéries sont Gram-positives ou Gram-négatives et la séquence d'acide nucléique est la sonde 1599, 1656 ou 1746.

FIGURE 1: THE THREE KINGDOMS
EUKARYOTES

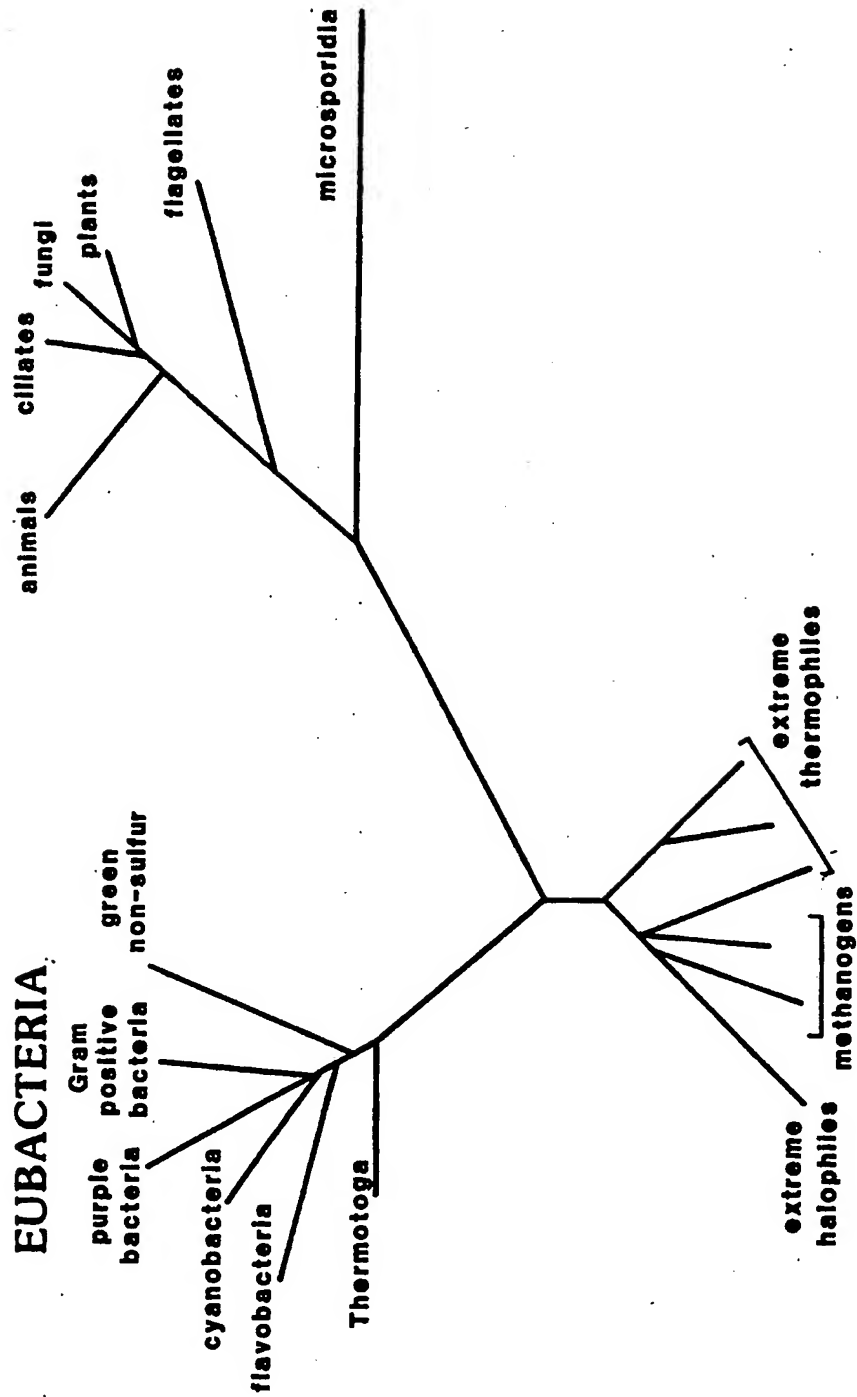


FIGURE2: THE EUBACTERIAL KINGDOM

